



*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C

## GC and GC/MS Q Exactive and applications

**Miloš Korman**


EMEA Sales Manager GC/GCMS

Sofia, October 22<sup>nd</sup>, 2016

The world leader in serving science

- **GC&GCMS portfolio 2016**
  - Expanded GC/GCMS product offering
- **GC news for 2016**
  - Injectors & Detectors
- **GCMS ISQ news for 2016**
  - GC/GCMS Product Portfolio and Application
- **GCMSMS TSQ news for 2016**
  - Application in Tox & Forensic
- **QEGC Update**
- **Conclusion**

# GC & GC-MS Products portfolio

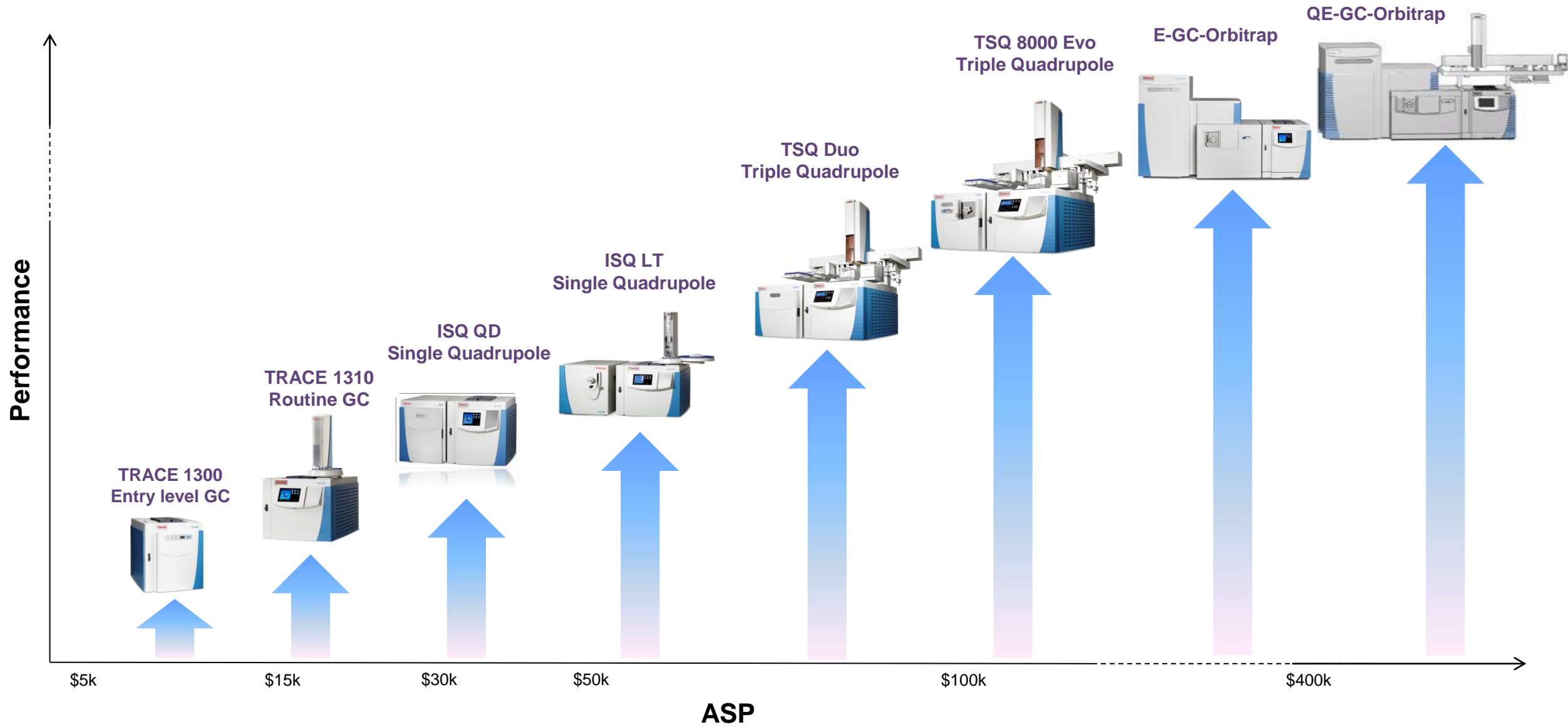



*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C

**GC news for 2016**

# TRACE 1300 Series GC – Two Models Available

**TRACE 1310 GC:** Touch screen interface provides instant access for ease of use and local control



**TRACE 1300 GC:** Local built-in ultra-simplified user interface – two buttons and four LEDs



## Modules available:

**Injectors:** SSL - SSL backflush - PTV - PTV backflush - **TSI** - **OCI** - **Helium Saver SSL+**

**Detectors:** FID - TCD - ECD - NPD – **FPD** - **PDID** - MS\*

**Other Options:** **MS NoVent**, **Auxiliary Oven**, Inj/Oven Cryo, Aux carrier

**Software drivers:** Chromeleon **6.8**, 7.1 and **7.2** CDS, Xcalibur, ChromQuest, ChromCard  
Agilent™ **OpenLAB™** EZChrom ed. CDS, Agilent **EZChrom™** (Agilent is a trademark of Agilent Corporation)

+ **Newly Released Modules and Options**

\*Supported Thermo Scientific MS: ISQ™ Series GC-MS; ITQ™ GC-MS; TSQ Quantum XLS GC-MS/MS; DFS; Delta V




*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C


**Injectors update for 2016**

# What is iC-TSI= Instant Connect ThermoSpray Injector

temperature profile comparison

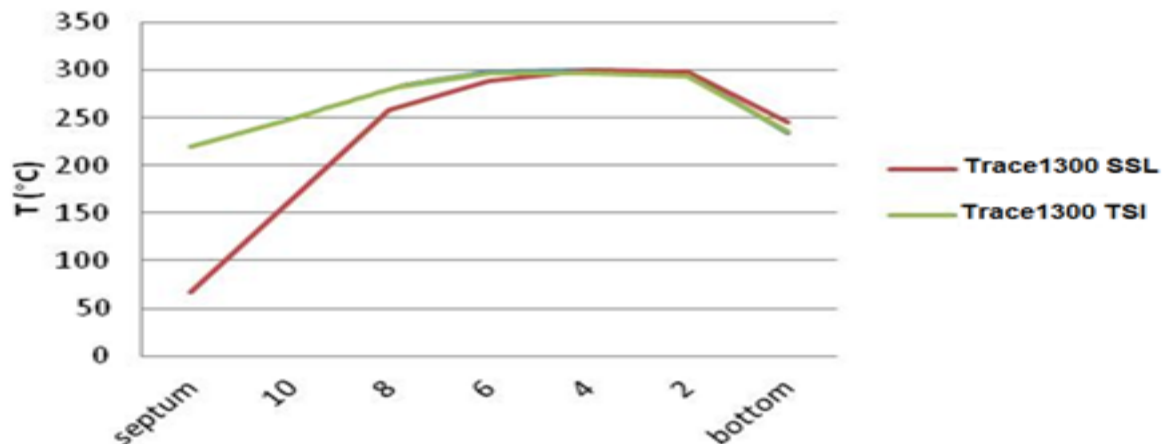

recovery comparison

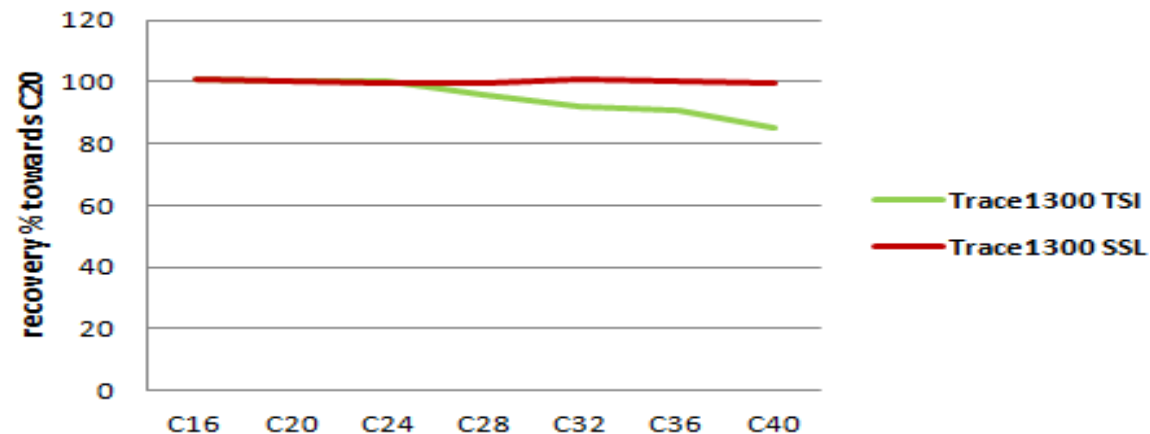

## TRACE1300 SSL

**Cold head**

**Liquid band injection**

**Needs fast injection**

**Simple and easy to set up**

**Not optimized for manual injection**

**Liner with glass wool required**

**Agilent like SSL injector**

## TRACE1300 TSI

**Hot head**

**Thermospray injection**

**Suitable for “Hot Needle Injection” technique**


**Method transfer from TRACE Ultra GC**

**Well performing with manual injection**

**Empty liner is used**

**Required for specific applications**

# i-TSI module

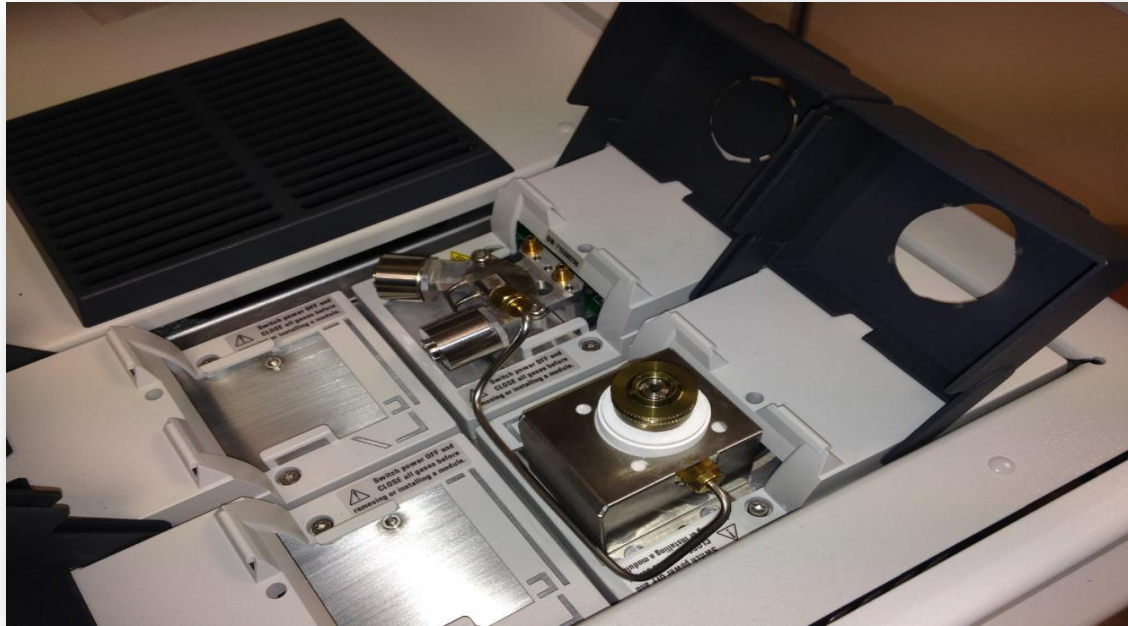

**p/n 19070016 TRACE 1300/1310 TSI MODULE**

It occupies two injector slots onto the Trace1300 GC Series



# Some of the challenging pesticides...

## ✓ Azole family (aim to tail in injector/T-line)

- Azaconazole
- Bitertanol
- Bromuconazole

## ✓ Pyrethrinoids (always difficult)

- Deltamethrine
- Cyhalothrines

## ✓ OrganoChlorine compounds

- Endosulfan sulfate

## ✓ Uracyl herbicide

- Bromacil (difficult detection)

## ✓ Organophosphorus compounds (sensitive to injector activity)

- Dichlorvos
- Vamidothion
- Demethon S methyl Sulfone
- Dimethoate

## ✓ Phenyl sulfamide insecticides

- Tolyfluanide


## ✓ Benzamide herbicides family

- Isoxaben (sensitive to column, T-line...)


# Isoxaben

i-SSL

Signal @10µg/L

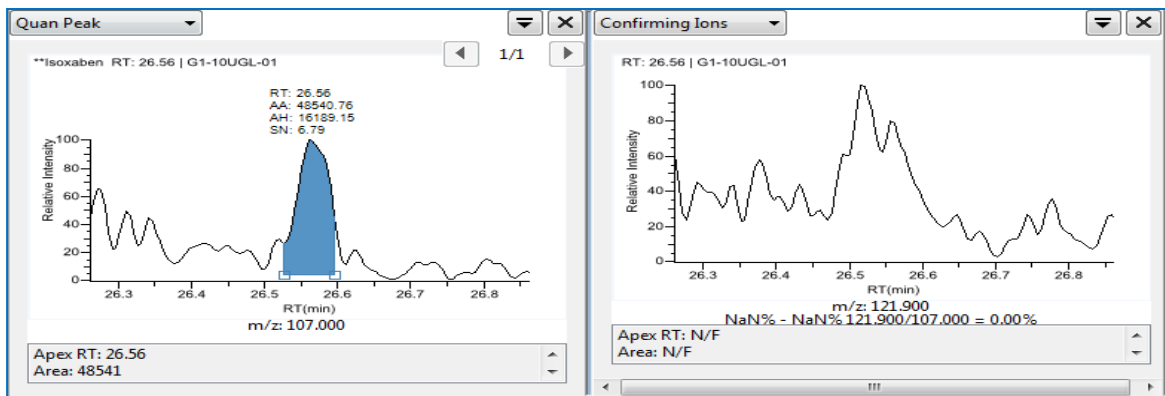


Linearity from 10pg to 30pg

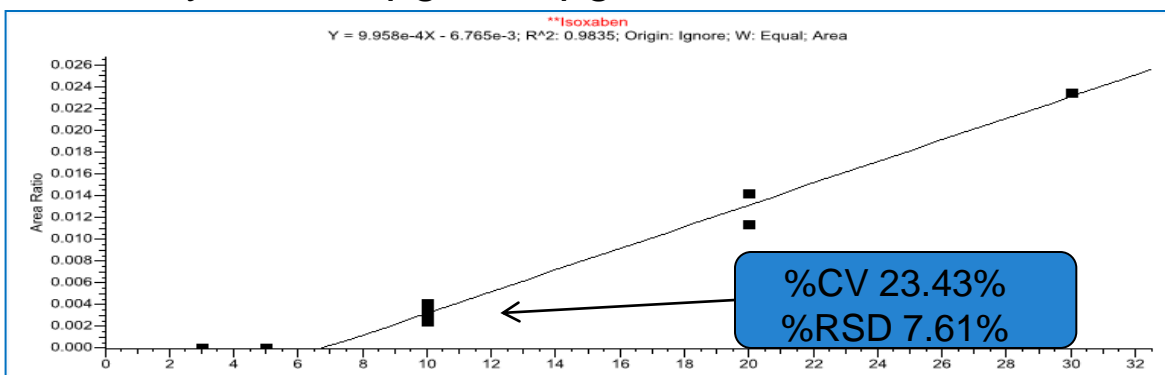


5 replicates @10µg/L

i-TSI

Signal @3µg/L

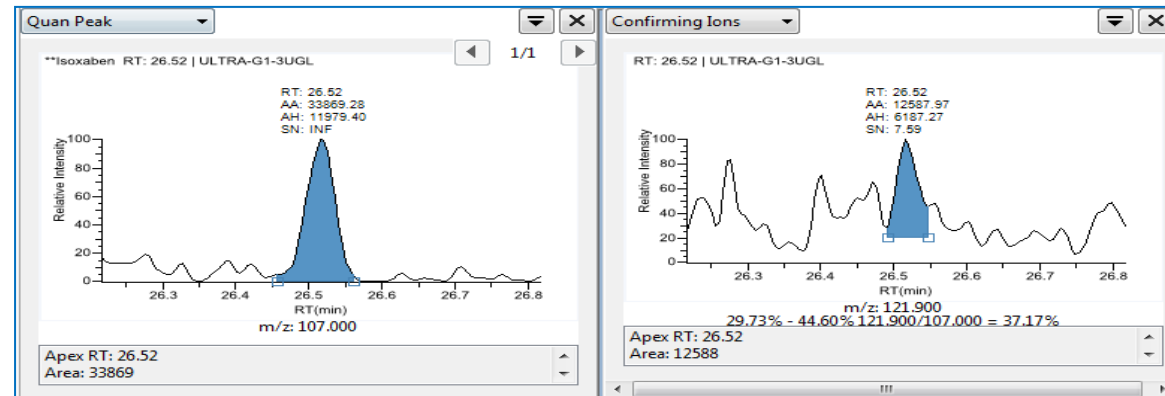

Linearity from 3pg to 30pg

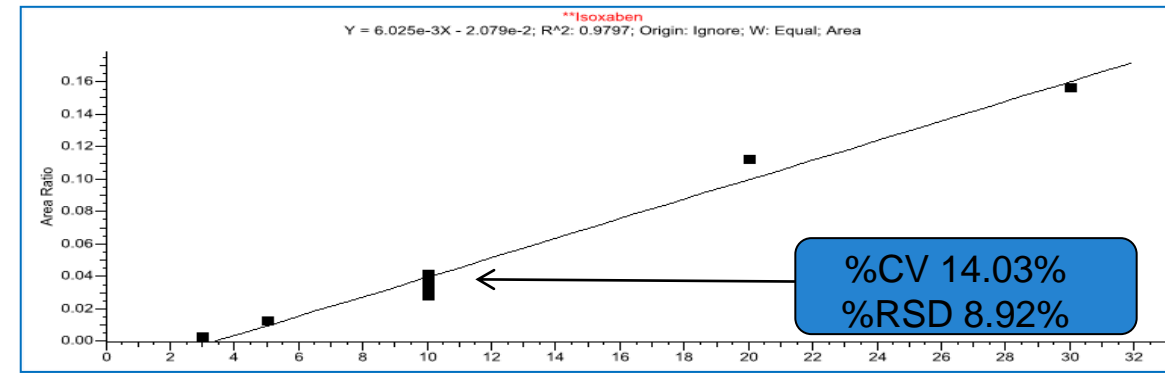


5 replicates @10µg/L


# Bromacil

i-SSL

Signal @10µg/L

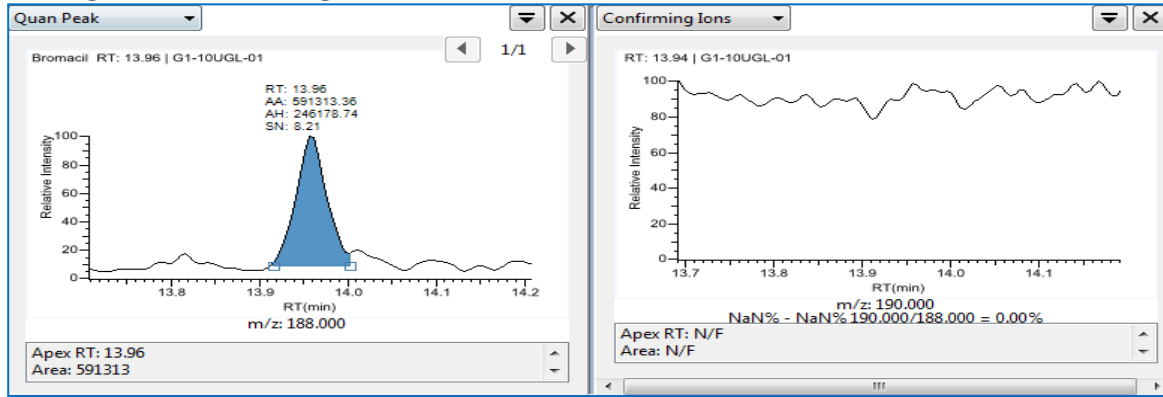


Linearity from 3pg to 30pg


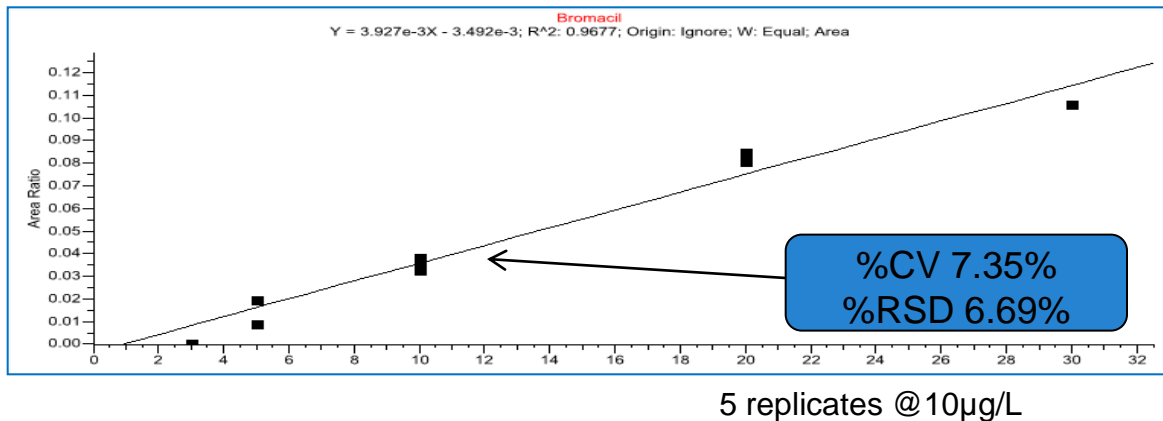


Signal @3µg/L

i-TSI

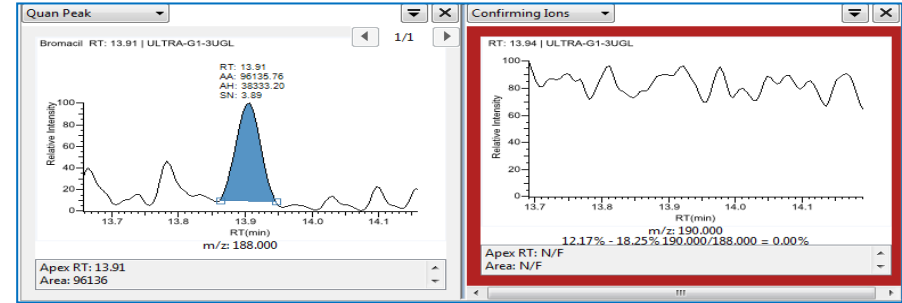


Signal @10µg/L

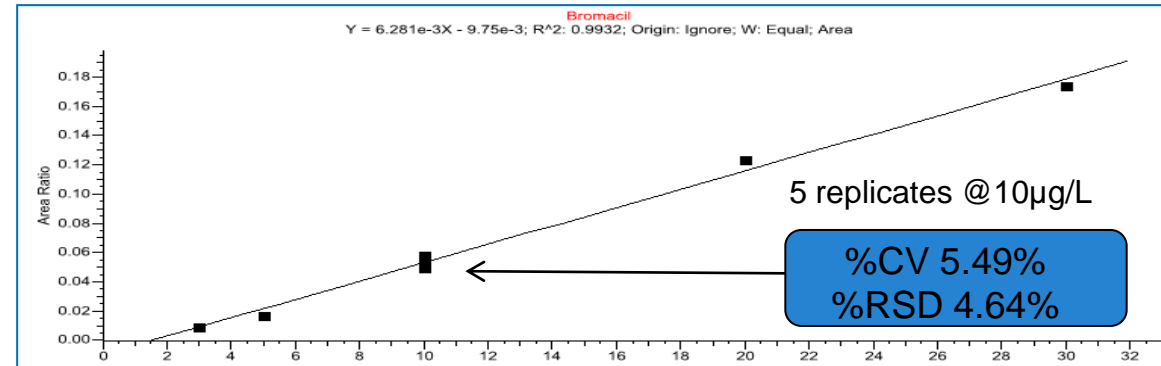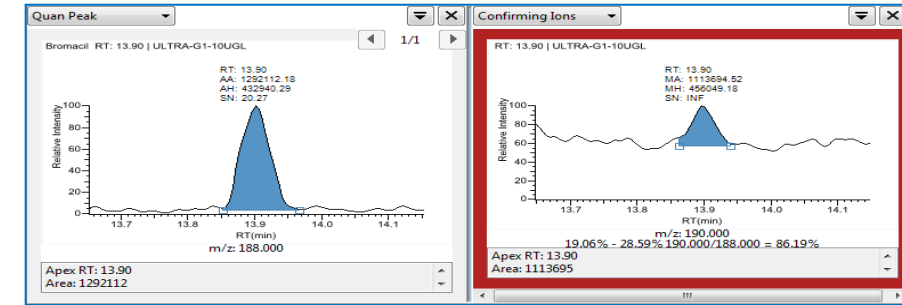

# Endosulfan sulfate


i-SSL

Signal @3 $\mu$ g/L

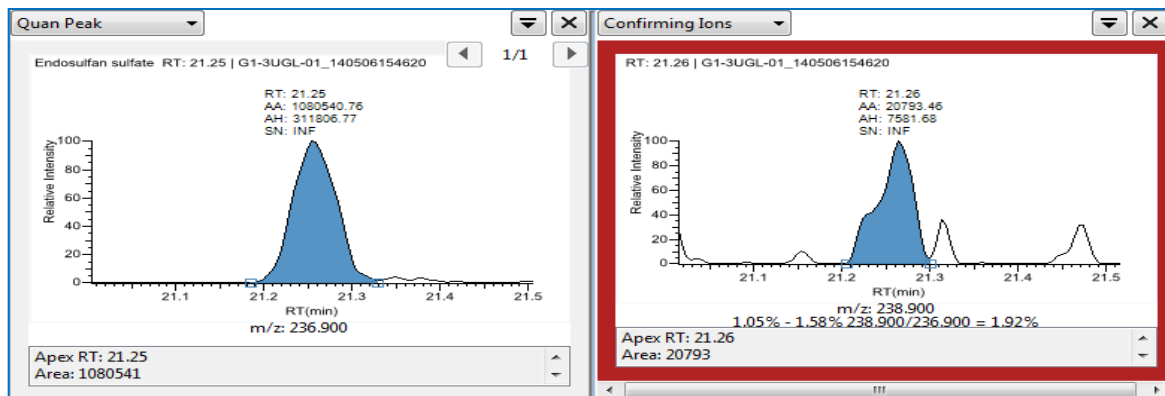


i-TSI


Signal @3 $\mu$ g/L

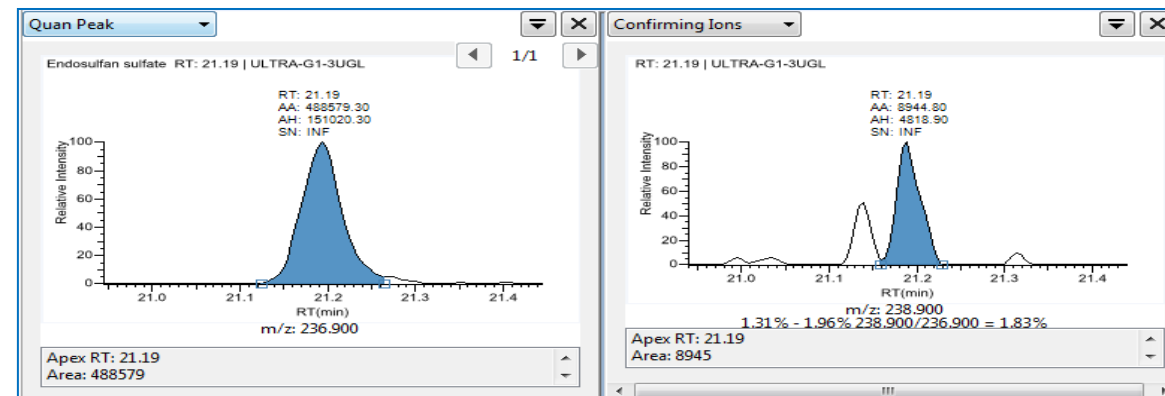

Linearity from 3pg to 30pg

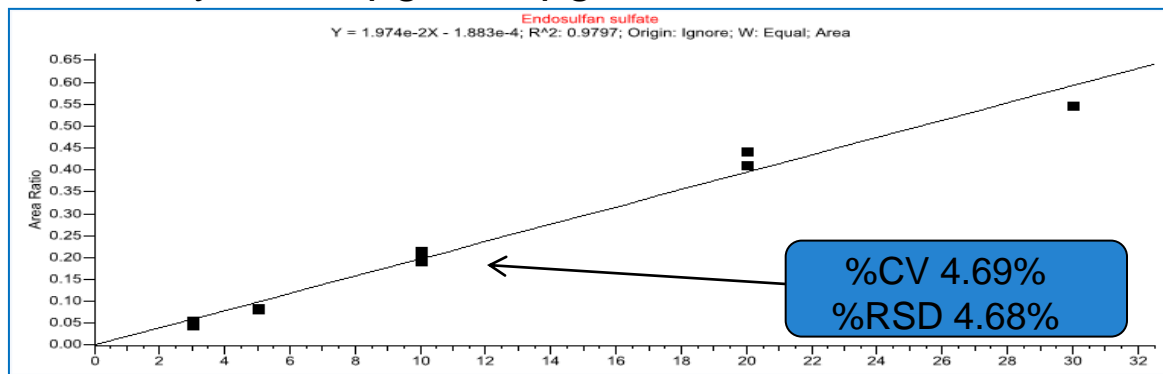

5 replicates @10 $\mu$ g/L

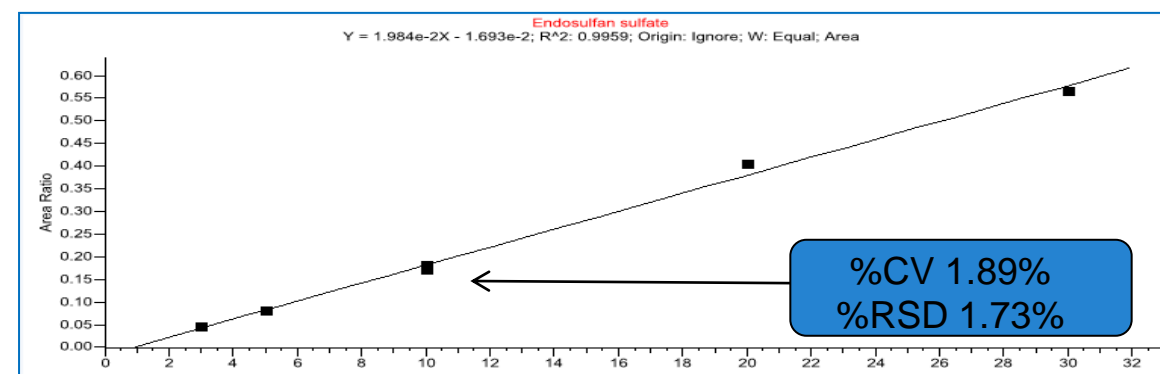

5 replicates @10 $\mu$ g/L

# When i-TSI is the right solution?

Looking at limit of detection, SN, area, %CV at 10pg and linearity we can observe that :

- ✓ i-TIS is a clearly better with improved limit of detection
  - ✓ Isoxaben
  - ✓ Vamidothion
  - ✓ Bromacil
  
- ✓ i-TIS is slightly better with higher SN
  - ✓ Endosulfan sulfate
  - ✓ Bitertanol
  - ✓ Cyfluthrines
  - ✓ Deltamethrines
  - ✓ tolyfluanid
  
- ✓ i-SSL and i-TIS are as excellent in matter of other Pesticides, detection, linearity and repeatability
  
- ✓ Thermospray injection can be advantageous with some critical compounds particularly prone to thermal degradation

# iC-OCI: Instant Connect On-column Injector module



iC-OCI inlet



- ✓ Based on the Carlo Erba cold on-column injector
- ✓ Merlin Microseal septum for low maintenance
- ✓ Continuous cooling of the injector body thanks to the integrated fan
- ✓ Pre-column secondary cooling
- ✓ Compatible with manual injection and Triplus RSH autosampler
- ✓ Two OCI modules can be mounted by using 1 set of Aux Modules

# iC-OCI: Instant Connect On-column Injector module

How to order:



Auxiliary Gas Module

Aux Temperature/Cryo Module



*iC-OCI Module kit*

<b>p/n 19071065</b>	<b>iC-OCI Module kit</b>	<i>Includes the Merlin Microseal 26s, the solenoid valve for secondary cooling with mounting bracket, tubing</i>
<b>p/n 19070080</b>	<b>Aux Carrier Gas Module</b>	<i>For Carrier Gas supply and control</i>
<b>p/n 19070070</b>	<b>Aux Temp Module</b>	<i>For the timed actuating of the secondary cooling</i>


# iC-OCI: the right solution

- ✓ Thermolabile compounds → PBDE 209
- ✓ Very high boilers → triglycerides, HT SimDist
- ✓ Wide boiling point range
- ✓ When an alternative to the PTV injector is required
- ✓ When sample does not include no-elutable by-products (clean sample extract)

## Benefits

- ✓ True cold injection → no thermal breakdown
- ✓ Direct injection into the column → no discrimination
- ✓ Sensitivity → splitless and LV injection
- ✓ Excellent injection repeatability < 0.5% RSD

Inj.	C12	C14	C16
1	1729812	1718526	1700241
2	1736398	1721253	1703718
3	1730131	1717829	1699899
4	1741246	1726306	1708989
5	1744551	1728775	1710122
6	1746308	1731325	1713732
7	1751230	1731861	1717433
8	1740056	1725927	1707828
9	1750468	1735191	1717330
10	1744906	1729495	1709880
Average	1741510	1726649	1708917
Stdev	7551	5851	6269
RSD %	0,43%	0,34%	0,37%

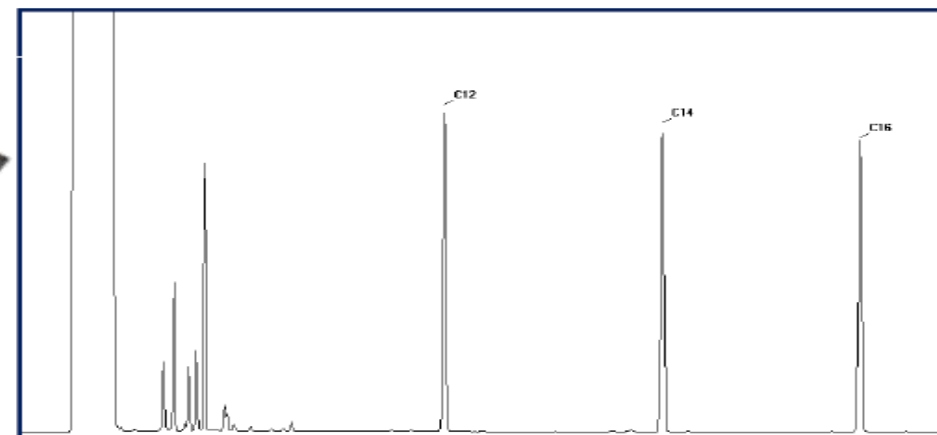




## Triglycerides

Blue = Blank hexane


Pink = Diet-type margarine

Black = inesterified triacylglyceride mixture

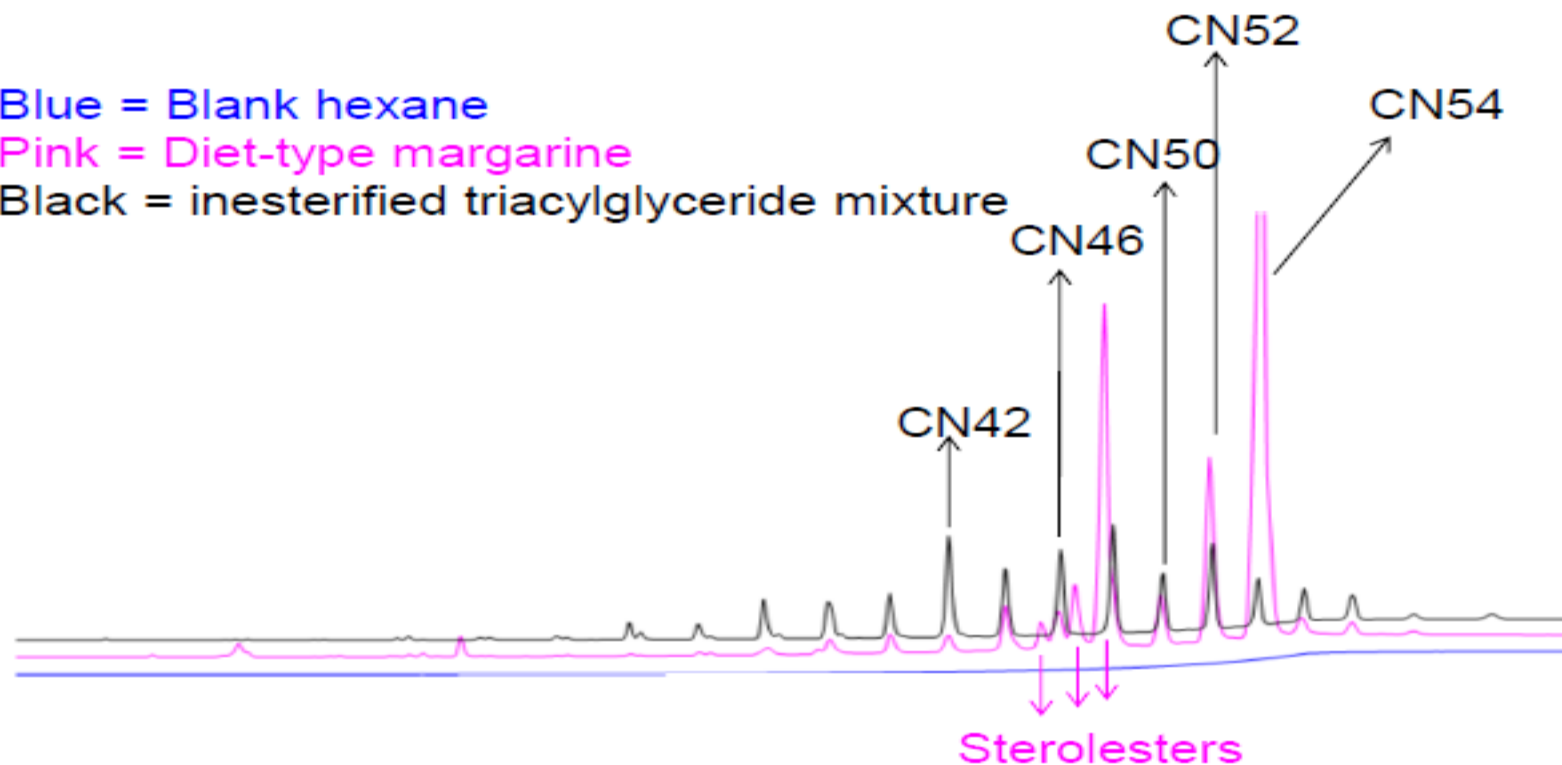

*CN54 = Tri-C18, glycerol backbone with three C18 fatty acids attached.*

## HT SimDist

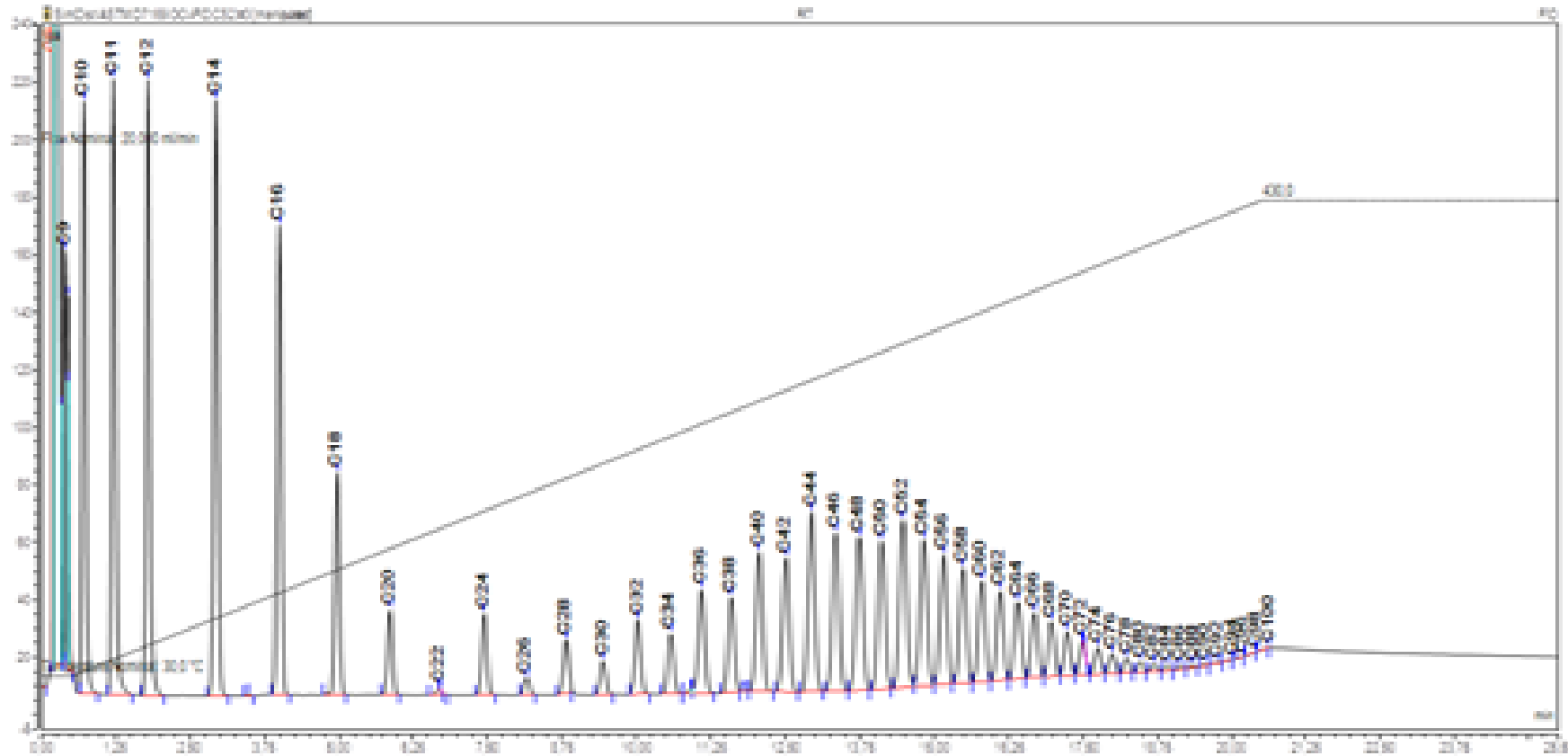



*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C

**Detectors update for 2016**


# Instant Connect FPD detector

- Available as Instant Connect module, can be mounted on any existing TRACE 1300 Series GC
- Highest sensitivity for the determination of Sulfur, Phosphorous and Tin containing species
- Acquisition rate up to 300 Hz
- Available, as option, as dual wavelength option with two FPD independent Instant Connect modules



# iC-PDD Pulse discharge detector

- Alternative to universal detector for ppb concentration range
- Engineered solutions (eg. specialty gases)
- Alternative for ECD if radioactive source is not suitable
- Replacement for FID, PID specifically required for the low detectability range

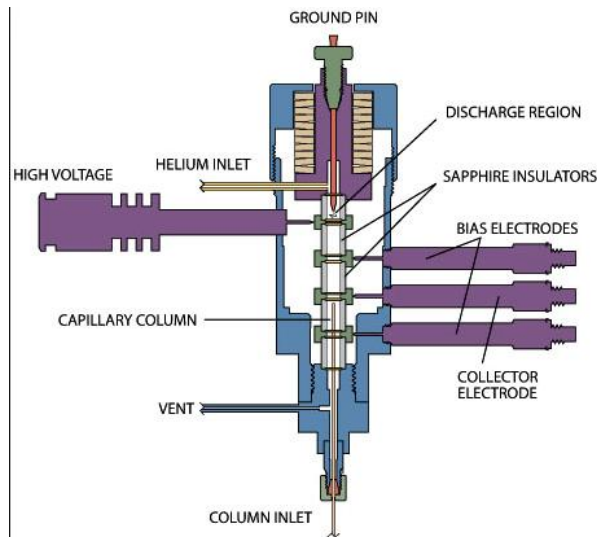

Instant Connect Pulsed Discharge Detector (PDD) for TRACE™ 1300 GC Series

# PDD minimum detectable amounts

MDQ S/N=3

*Sample injected 250uL, 1/10 split ratio*

<b>compound</b>	<b>pg</b>	<b>ppb (v/v)</b>
Carbon dioxide	0,9	30
Carbon monoxide	4	150
Nitrous oxide	0,9	30
Methane	0,6	20
Acetylene	1	40
Ethylene	0,7	25
Ethane	0,7	25






*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C

**Autosamplers update for 2016**

# Thermo GC Autosamplers Lineup

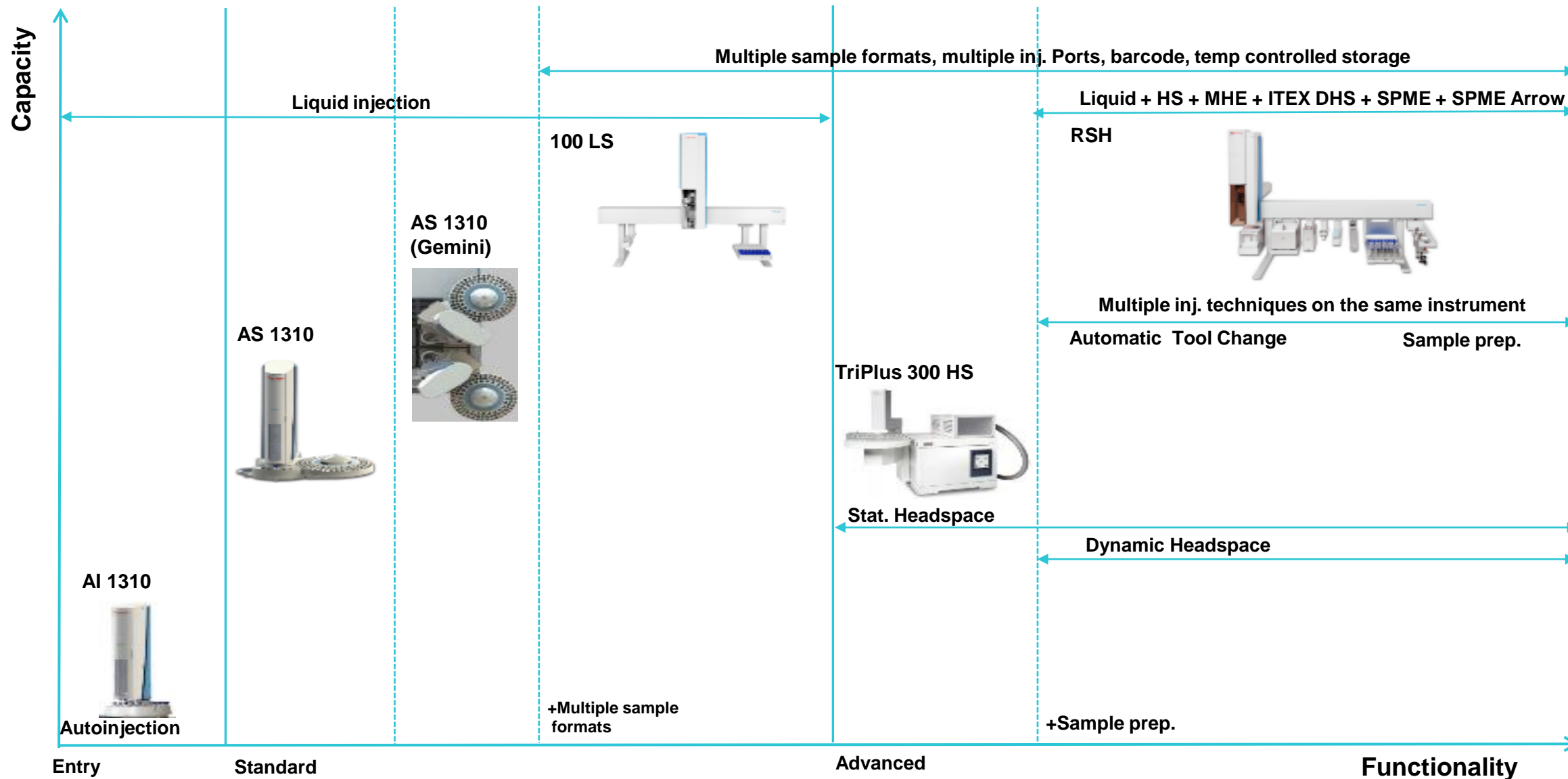



# Solvent-free extraction of volatiles


## *No enrichment of volatiles*

- Static Headspace Sampling (SHS)

## *With enrichment of volatiles on adsorbent material*


- Purge&Trap (P&T)
- Solid Phase Micro Extraction (SPME, MEPS)
- Stir Bar Sorptive Extraction (SBSE)
- Headspace Sorptive Extraction (HSSE)
- Dynamic Headspace Sampling (DHS, ITEX)
- Thermal Desorber (TD)

# Triplus RSH expanded options for VOC and SVOC market

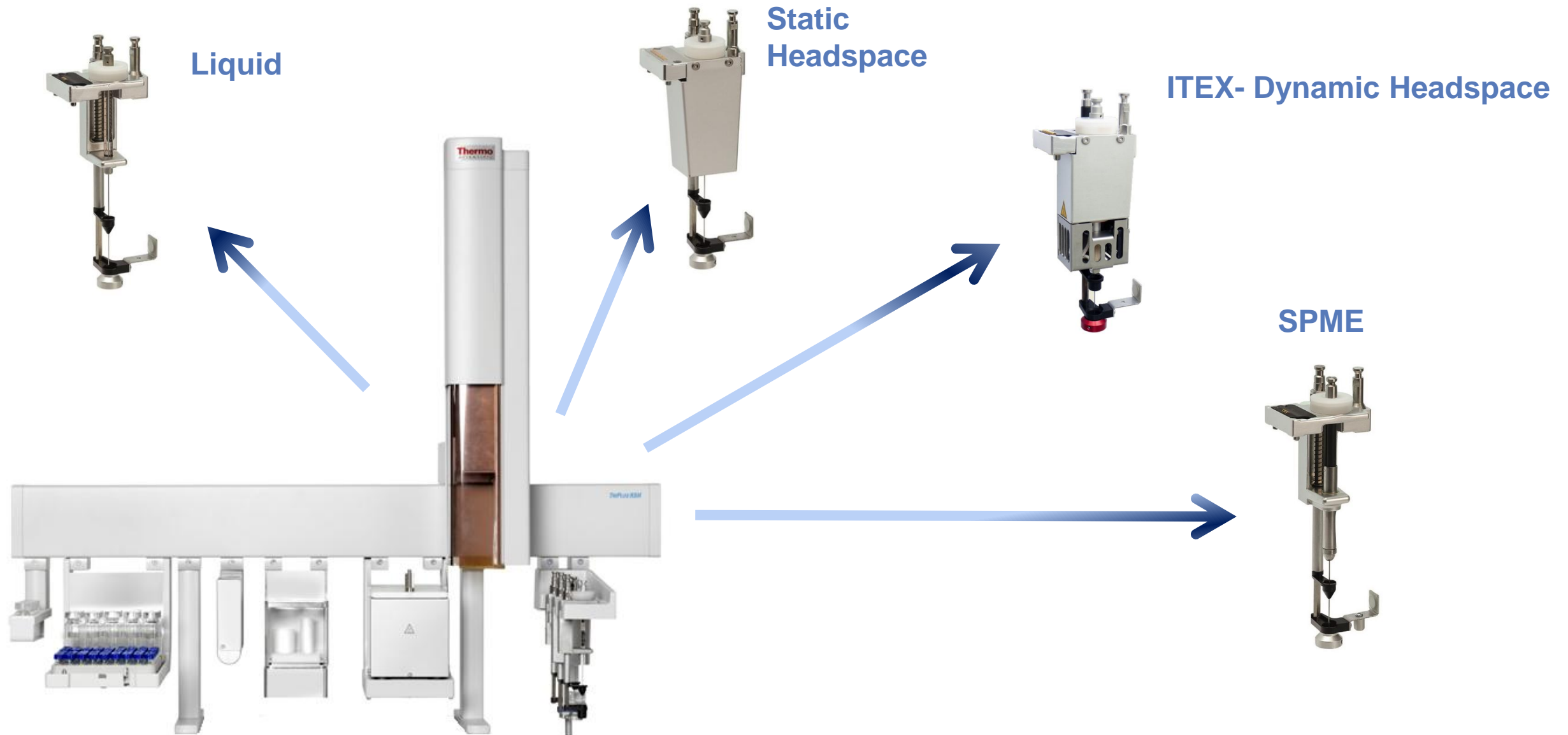


# ITEX Dynamic Headspace

- ITEX Workflow

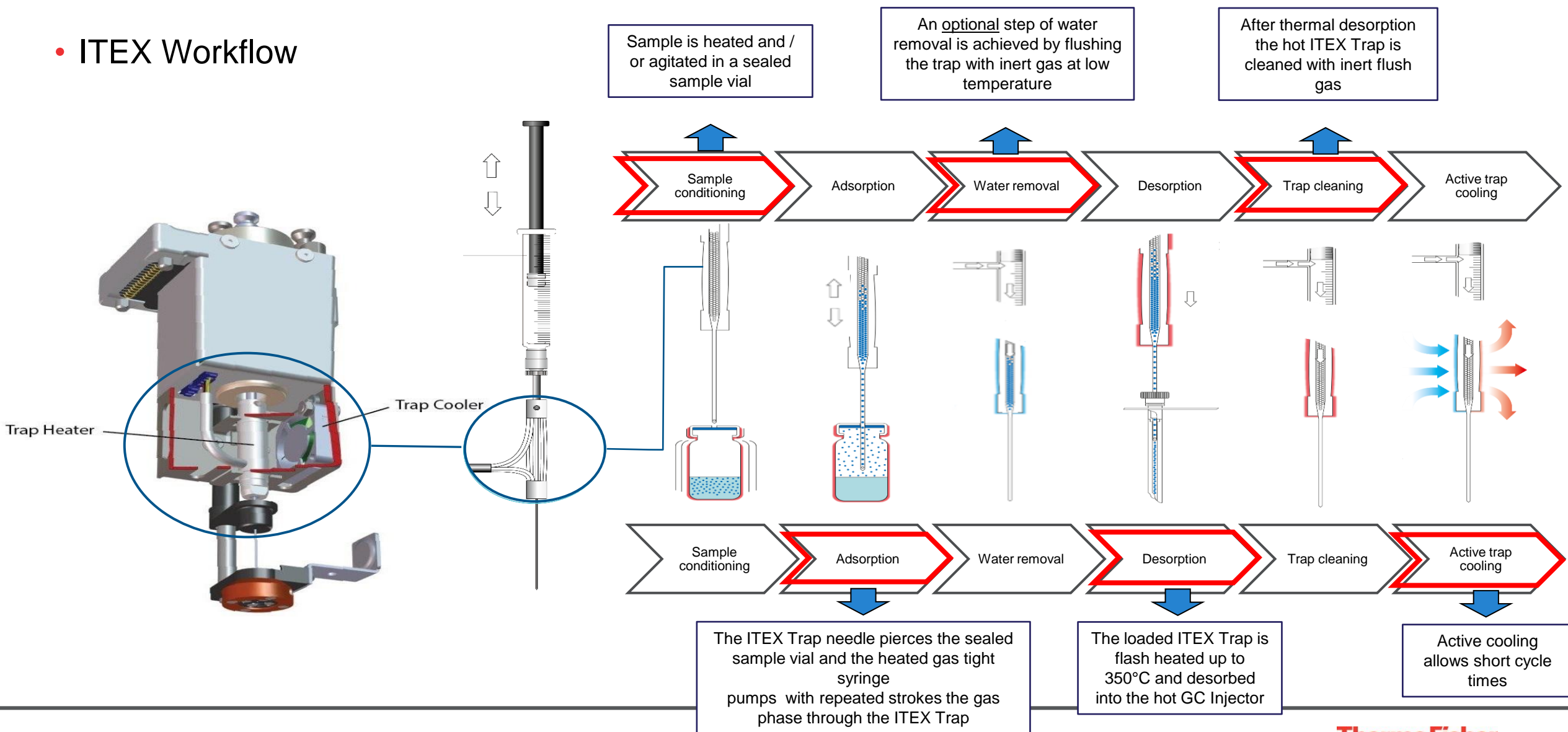


## .....Simplicity and handiness

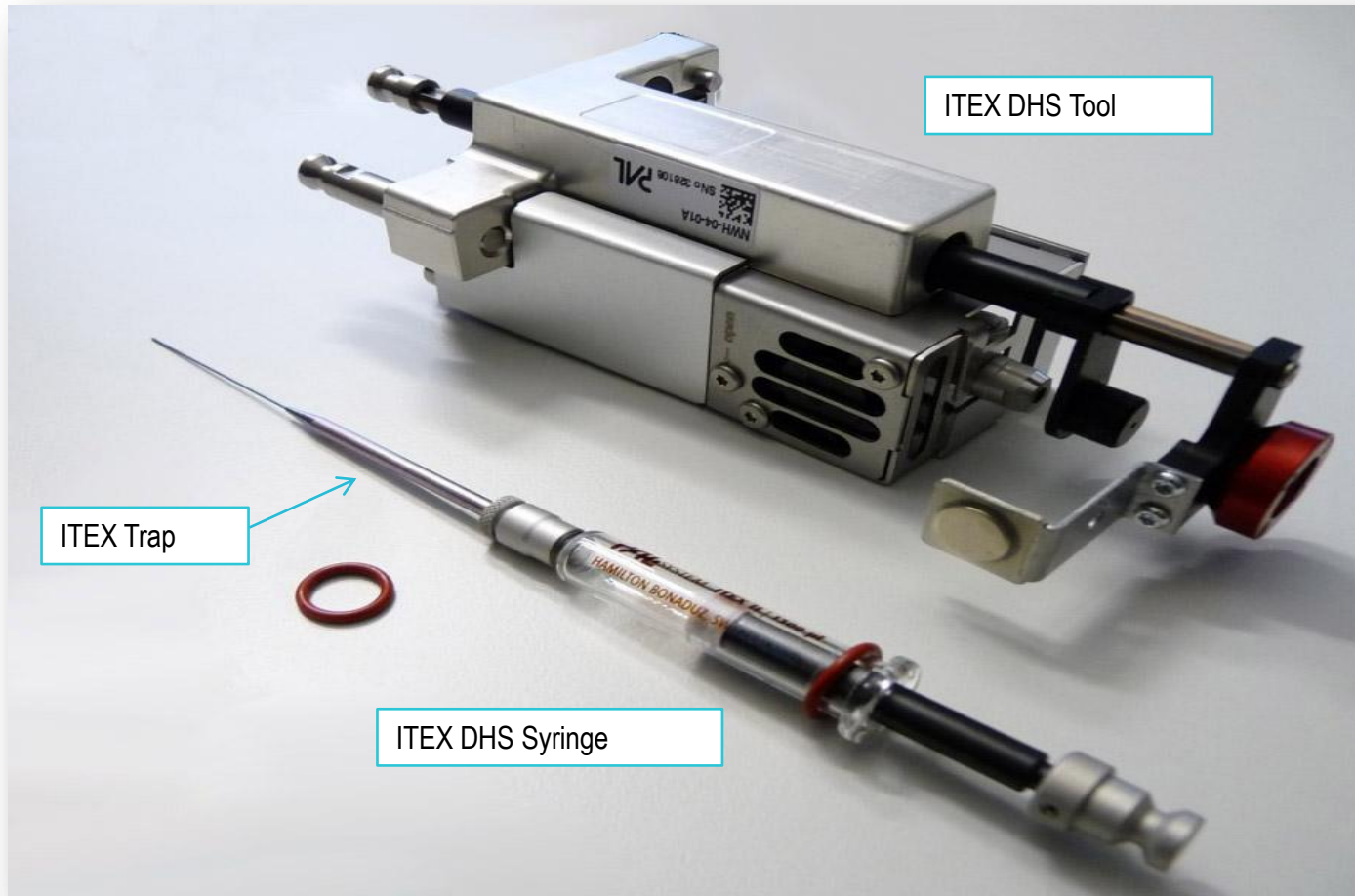

- Easy sample handling
- No loop, transfer line, or switching valve
- Easy access to the trap and syringe
- Straightforward troubleshooting and maintenance
- Minimum instrument downtime

## Proven industry standard adsorbents

- Standard ITEX DHS adsorbent:
  - Tenax TA 80/100 mesh (volatile and semivolatile compounds)
- Additional ITEX DHS adsorbents:
  - Carboxen C, 80/100 mesh
  - Carboxen 1000, 60/80 mesh
  - Carbosieve S III, 60/80 mesh
  - Tenax GR, 80/100 mesh
  - Tenax GR, 80/100 mesh / Carbosieve S III, 60/80 mesh

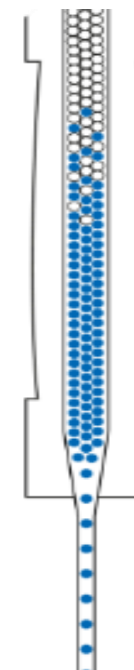


### **SUGGESTED APPLICATIONS:**

Tenax TA:	Volatile and semivolatile compounds, temperature limit of 350°C
Carboxen:	Hydrophobic properties minimized sample displacement by water
Carbosieve/Carboxen:	For very volatile compounds, e.g. Vinylchloride, Freon compounds

Reference: **Optimization strategies of in-tube extraction (ITEX) methods**  
*J. Laaks et al. Anal Bioanal Chem (2015) 4017:6827-6838*

# TriPlus RSH: ATC Automatic Tool Change


- Multiple injection methods on one instrument: Flexibility and Productivity
- Intelligent Switch to Dynamic Headspace if Static Headspace is not sensitive enough

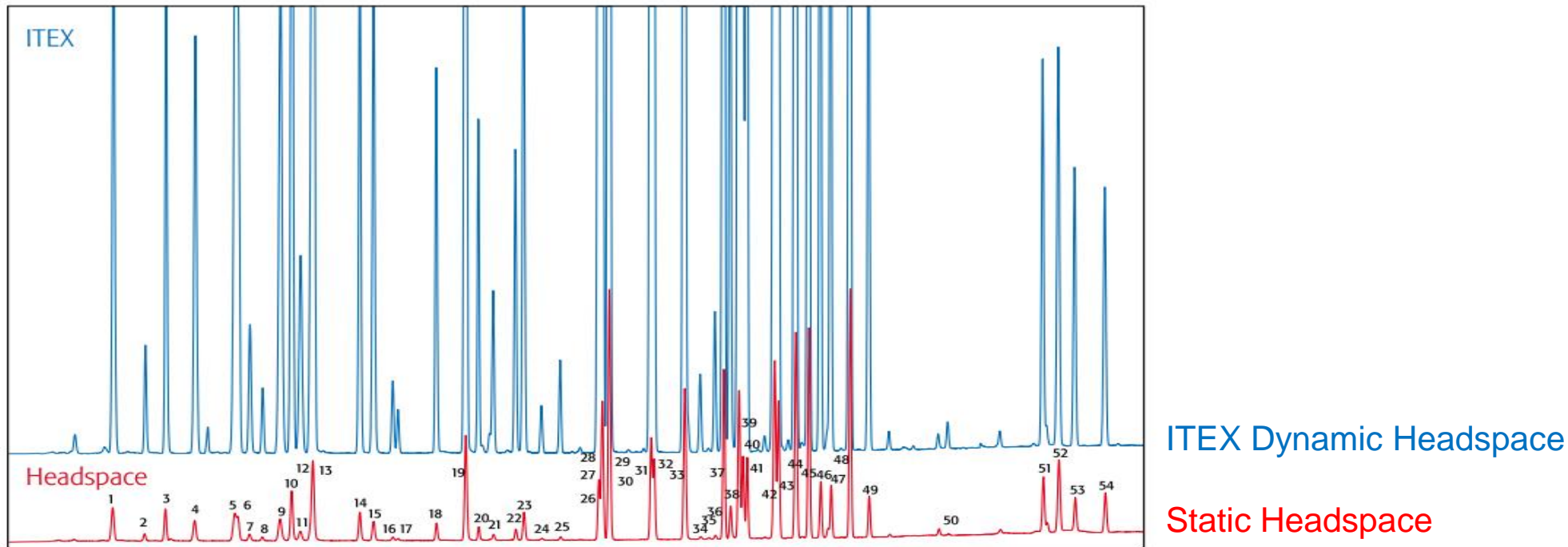

EPA 502.2: Increased sensitivity towards Headspace by magnitudes, **Easily reaching threshold values of US-EPA Methods for VOC's**

# TriPlus RSH: ATC Automatic Tool Change


- Switching of sample injection techniques on the same instrument
- TriPlus RSH is compatible with a wide range of syringe volumes
- ATC station enables **unattended** syringe exchange **during a sequence**

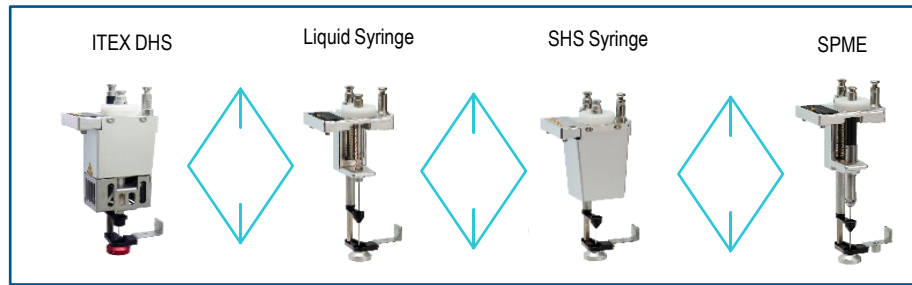


ATC (Automatic Tool Change)



TriPlus RSH: ATC (Automatic Tool Change)

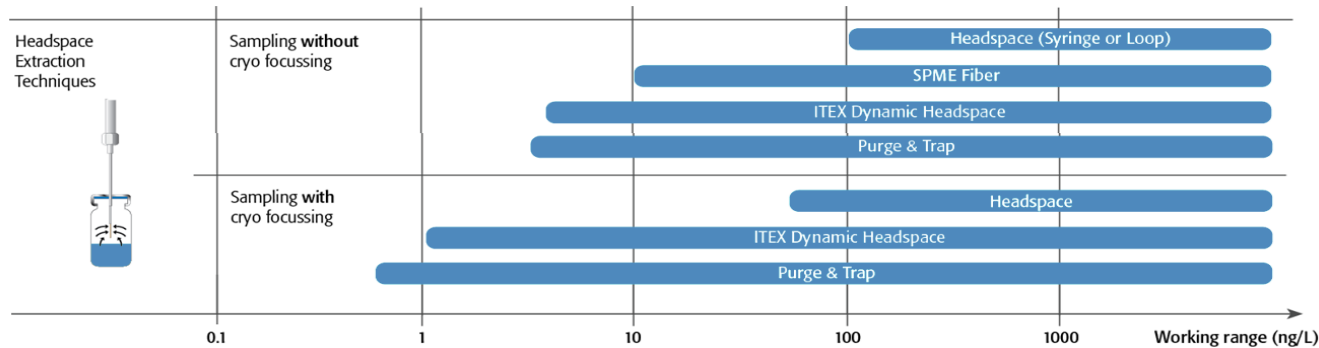


Achievable sensitivities by sampling technique

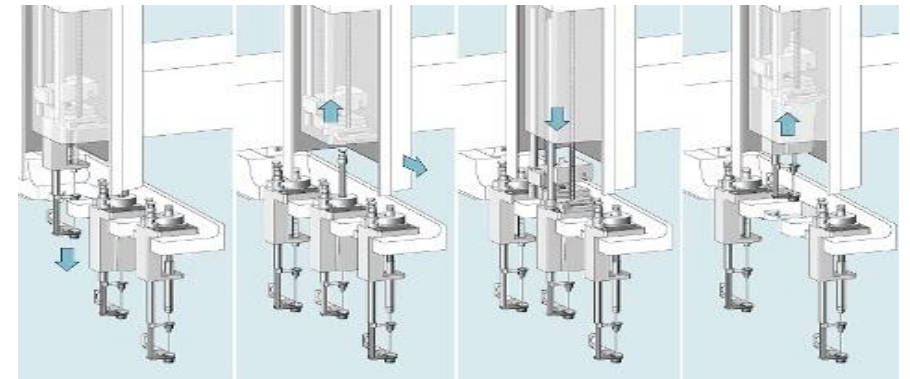


ATC (Automatic Tool Change)

# Triplus RSH

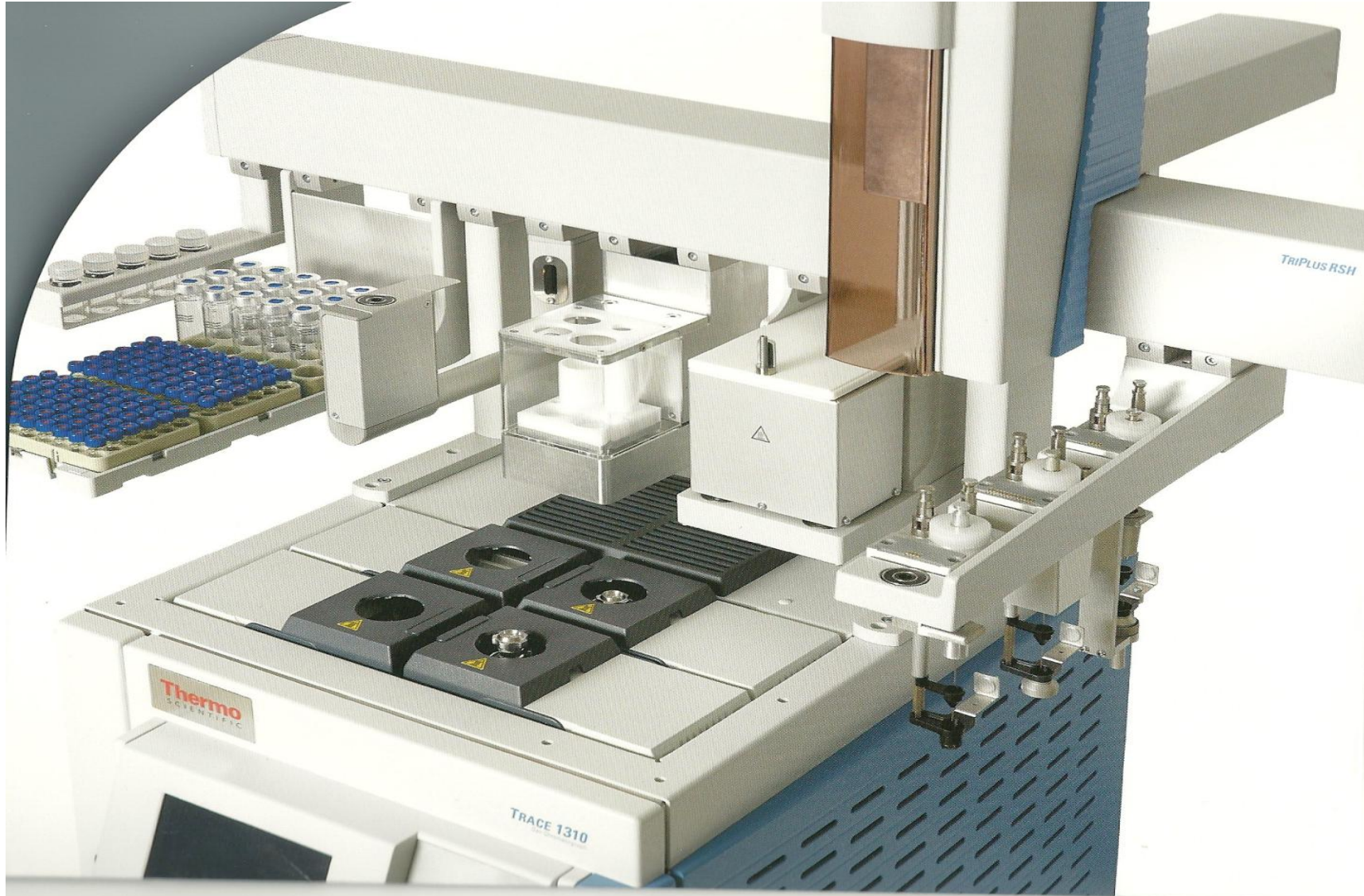



# Triplus RSH expanded options for VOC and SVOC market

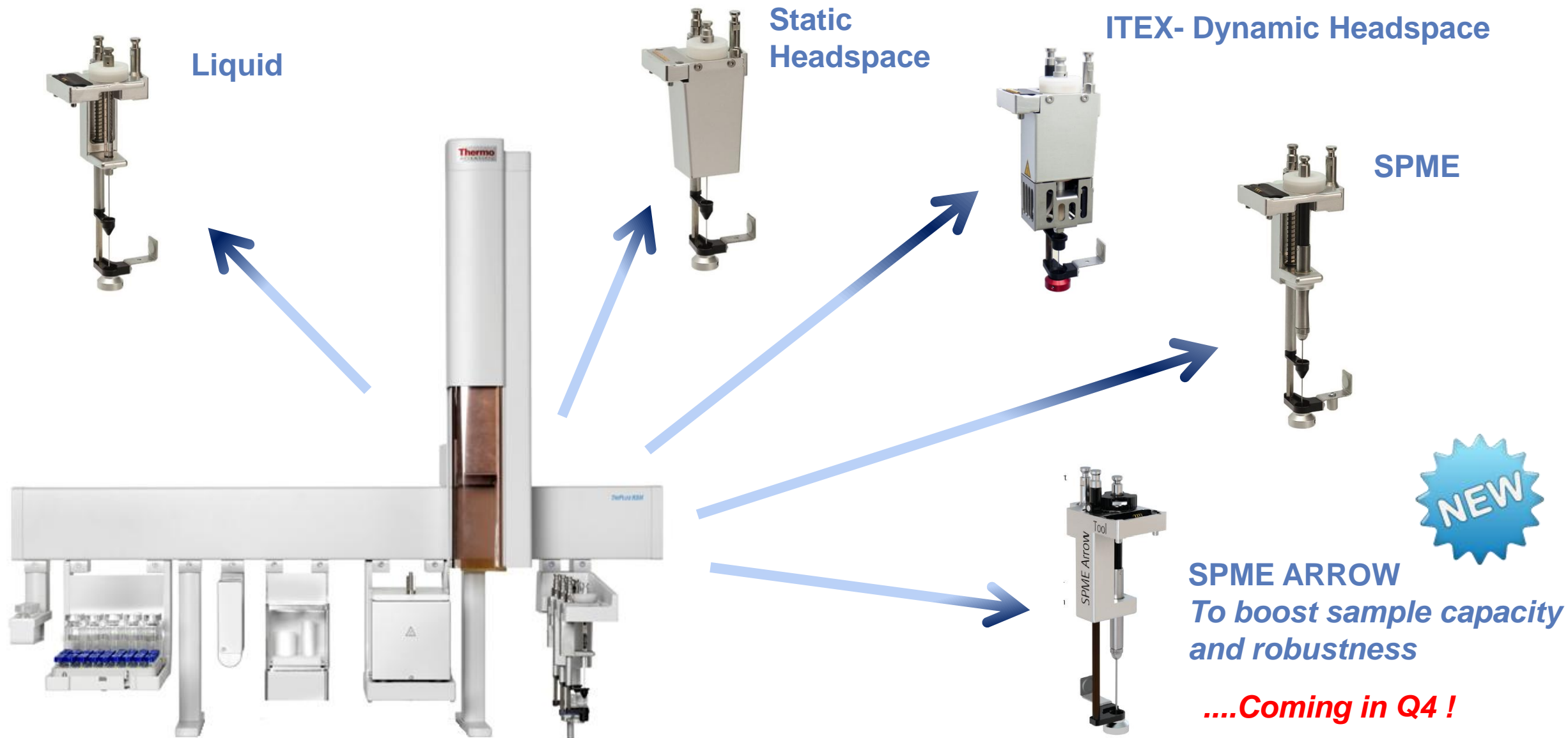

# SPME Arrow tools

- SPME Arrow tool




- Heatex-Stirrers



- SPME Arrow conditioning module




# SPME Arrow Fiber collection

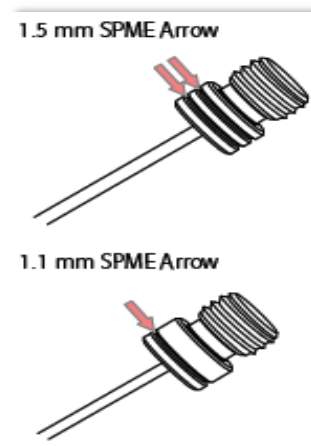

- A full range of fibers is available with different diameter and thickness
- Diameter is coded
- All SPME Arrow fibers have 20 mm phase length

Color code	Phase	Diameter	Phase Thickness
Red	PDMS	1,1 mm	100 um
Grey	Acrylate	1,1 mm	100 um
Light Blue	Carbon WR/PDMS	1,1 mm	120 um
Violet	DVB/PDMS	1,1 mm	120 um
Black	PDMS	1,5 mm	250 um

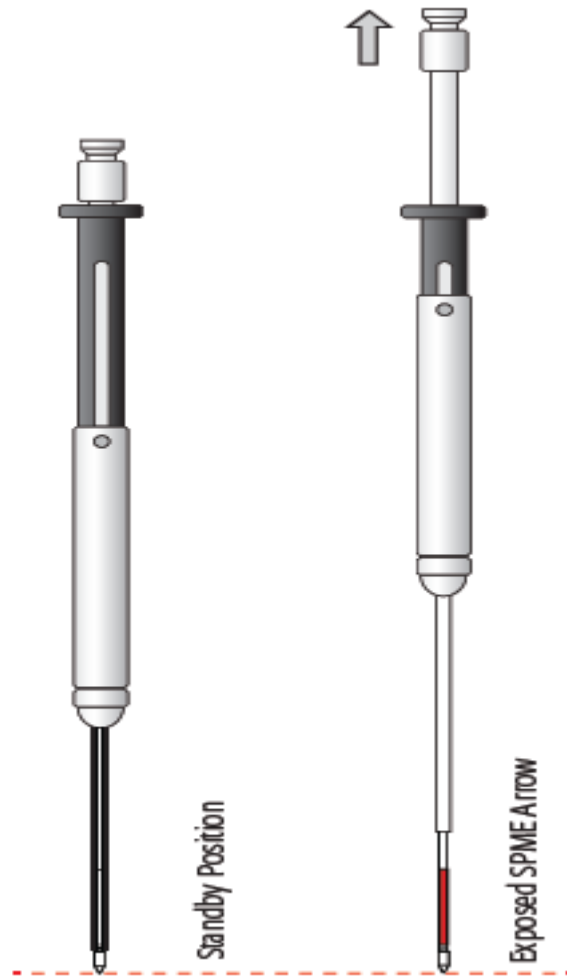


# Increased robustness



- Due to the larger needle diameter, an SPME Arrow adaptation kit for the Split/Splitless injector is necessary:
- Adaptation kit for TRACE 1300 Series GC Split/Splitless injector
- Adaptation kit for TRACE Ultra GC Split/splitless injector
- A specific liner (1.8 or 1.3 mm i.d. depending on the fiber diameter) is suggested for optimal performance


Higher volume = Up to 10 times more sensitivity

The table below shows the dimension of a PAL SPME Arrow 1.5 mm (a), 1.1 mm (b) and a SPME Fiber (c) in comparison:

		Sorption phase surface	Sorption phase volume
a		62.8 mm <sup>2</sup>	11.8 μL
b		44.0 mm <sup>2</sup>	3.8 μL
c		9.4 mm <sup>2</sup>	0.6 μL



# SPME Arrow

- Larger surface, faster extraction >> *2x higher sample throughput*
- More sorption phase, superior sensitivity>>*up to 10 times more sensitivity, wider dynamic range*
- Optimized geometry, greater robustness>>*SPME Arrow and injector septa last at least 2 x times longer*
- Designed for reliable automation via TriPlus RSH autosampler

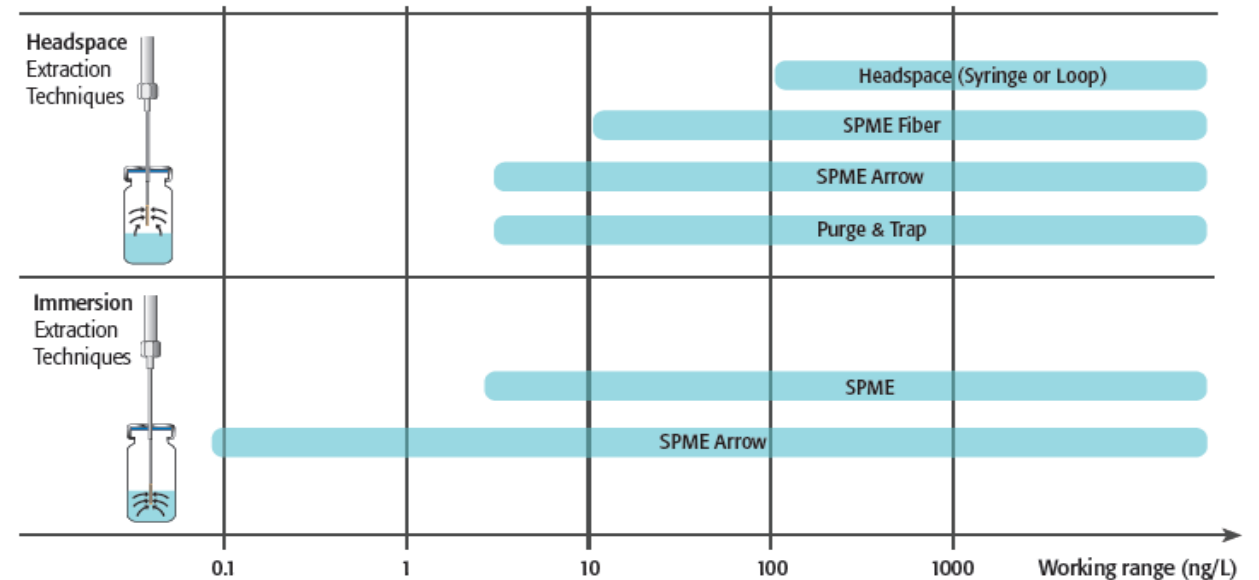

# SPME Arrow vs SPME fiber

- Headspace extraction: aroma analysis in white wine

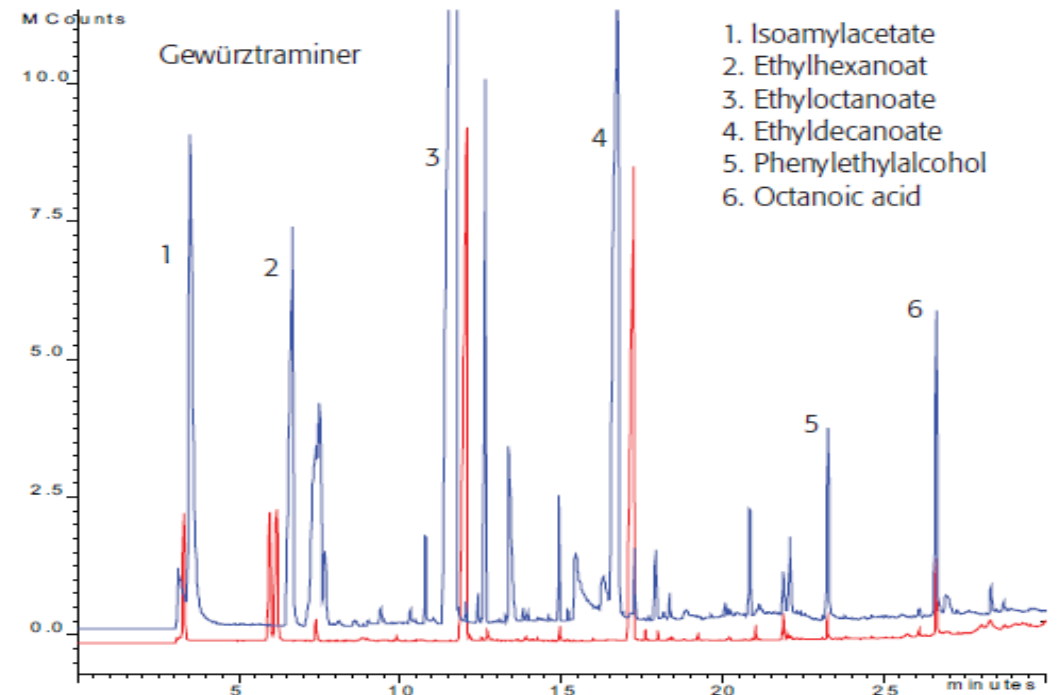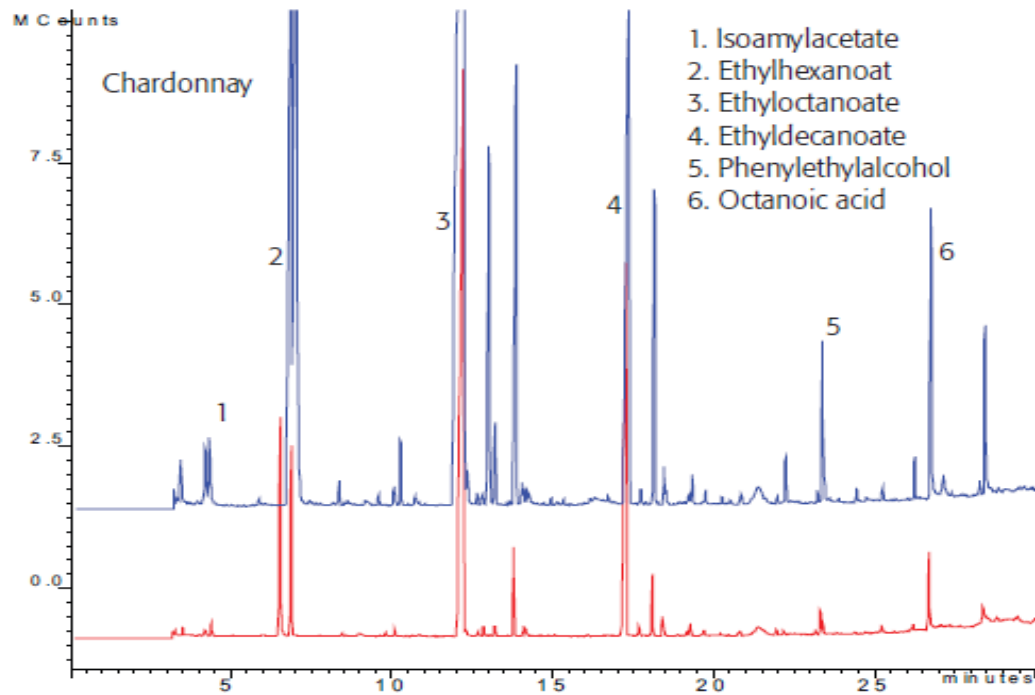

(SPME Arrow 100  $\mu\text{m}$ , 20 x 1.1 mm compared to SPME 100  $\mu\text{m}$ , 10 x 0.3 mm)

# SPME Arrow vs SPME vs SBSE

MDL and RSD obtained with SPME Arrow (250 umx 20 mm) for PAHs in water in comparison with literature data for classical SPME and SBSE bars

Compound	PAL SPME Arrow		SPME (Cheng et al.) [15]		SBSE (Carrera et al.) [22]	
	MDL (ng L <sup>-1</sup> )	RSD (%) (at 10 ng L <sup>-1</sup> )	LOD (SD X 3)	RSD (conc. at S/N=3×3)	LOD (conc. at S/N=3×3)	RSD (%) (at 50 ng L <sup>-1</sup> )
Naphthalene	0.3	5.7	2.7	9.0	/	/
Acenaphthylene	0.2	6.0	1.8	6.0	0.1	/
Acenaphthene	0.1	7.1	0.9	3.0	/	/
Fluorene	0.2	5.6	3	10.0	0.1	8.3
Phenanthrene	0.2	5.5	2.1	7.0	0.1	1.1
Anthracene	0.3	7.6	2.1	7.0	0.2	2.1
Pyrene	0.2	6.4	3.6	12.0	0.2	/
Fluoroanthene	0.2	6.2	2.1	7.0	0.2	/
1,2-Benzanthracene	0.1	6.2	2.1	7.0	0.2	6
Chrysene	0.1	11.0	1.5	5.0	0.2	10.6
Benzo(b)fluoroanthene	0.2	10.5	2.7	9.0	0.1	/
Benzo(k)fluoroanthene	0.2	8.6	1.8	6.0	0.1	/
Benzo[a]pyrene	0.3	7.2	3.6	12.0	0.1	/
Indeno(1,2,3 cd)pyrene	0.8	9.2	3.6	12.0	0.3	/
Dibenz(ah)anthracene	0.6	11.3	/	/	0.3	/
Benzo(ghi)perylene	0.8	11.9	1.8	6.0	0.3	/

MDL values calculated with a 99 % confidence interval


/ not determined

SPME Arrow - Evaluation of a Novel Solid-Phase Microextraction Device for Freely Dissolved PAHs in Water; Kremser A. et al 2016, 408, 943-952



# Immersion extraction: PAHs in water

Extraction yields for water samples containing 50 ng/L PAHs after 70 min

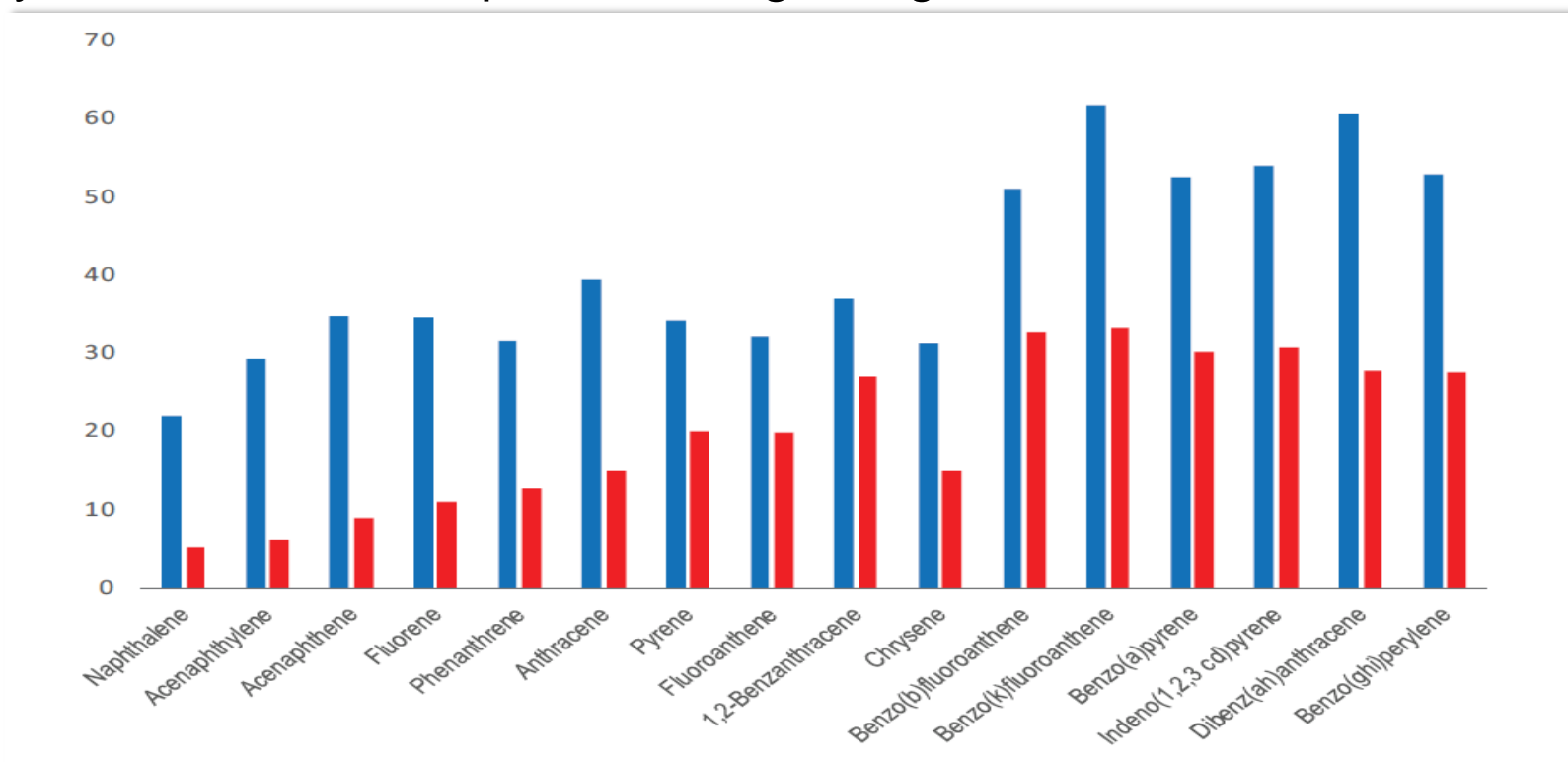

**SPME Arrow** 20 mm x 250 um carbon WR

**SPME fiber** 10 mm x 100 um Carbon WR

SPME Arrow - Evaluation of a Novel Solid-Phase Microextraction Device for Freely Dissolved PAHs in Water; Kremser A. et al 2016, 408, 943-952

## Sampling techniques comparison chart

Analytical quality	HS			DHS/Purge & Trap		SPME	SPME Arrow	ITEX-DHS	SBSE/HSSE
	Loop	syringe	P bal	headspace	sparkling	Immersion/ headspace	Immersion/ headspace	headspace	Immersion/ headspace
Sensitivity	+	+	+	++	+++	++	++++	+++	++++
Regulation compliance	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes
EPA methods	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	No
EU methods	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	No
Ease of use	+++	+++	+++	++	++	+++	+++	+++	+++
Sample handling	+++	+++	+++	++	+	+++	+++	+++	++
Analytes discrimination	low	low	low	low	med	med/high	med/high	low	high
Risk of contamination	low	very low	low	low	high	very low	very low	low	low

Hardware related	HS			DHS/Purge & Trap		SPME	SPME	ITEX-DHS	SBSE/HSSE
	Loop	syringe	P bal	headspace	sparkling	Immersion/ headspace	Immersion/ headspace	headspace	Immersion/ headspace
Modularity	No	Yes	No	No	No	Yes	Yes	Yes	Yes
Cost of device	high	low	high	high	high	low	low	low	high
Automated	Yes	Yes	Yes	limited	Yes	Yes	Yes	Yes	Yes
Space requirement	high	low	high	high	high	low	low	low	low
GC installation needs	high	low	high	high	high	low	low	low	low
GC injector blocked	Yes	No	Yes	Yes	Yes	No	No	No	No
Flexible w injection modes	No	Yes	No	No	No	Yes	Yes	Yes	Yes

# Head Space TriPlus 300 –transfer line




- **Instrument configuration**
  - TRACE 13XX GC
  - Dual column system
  - TriPlus 300 Head Space Autosampler
  - FID
  - Chromeleon 7.2 CDS



# Dual detectors confirmation analysis

- Screening with GC detector; confirmation with MS detector

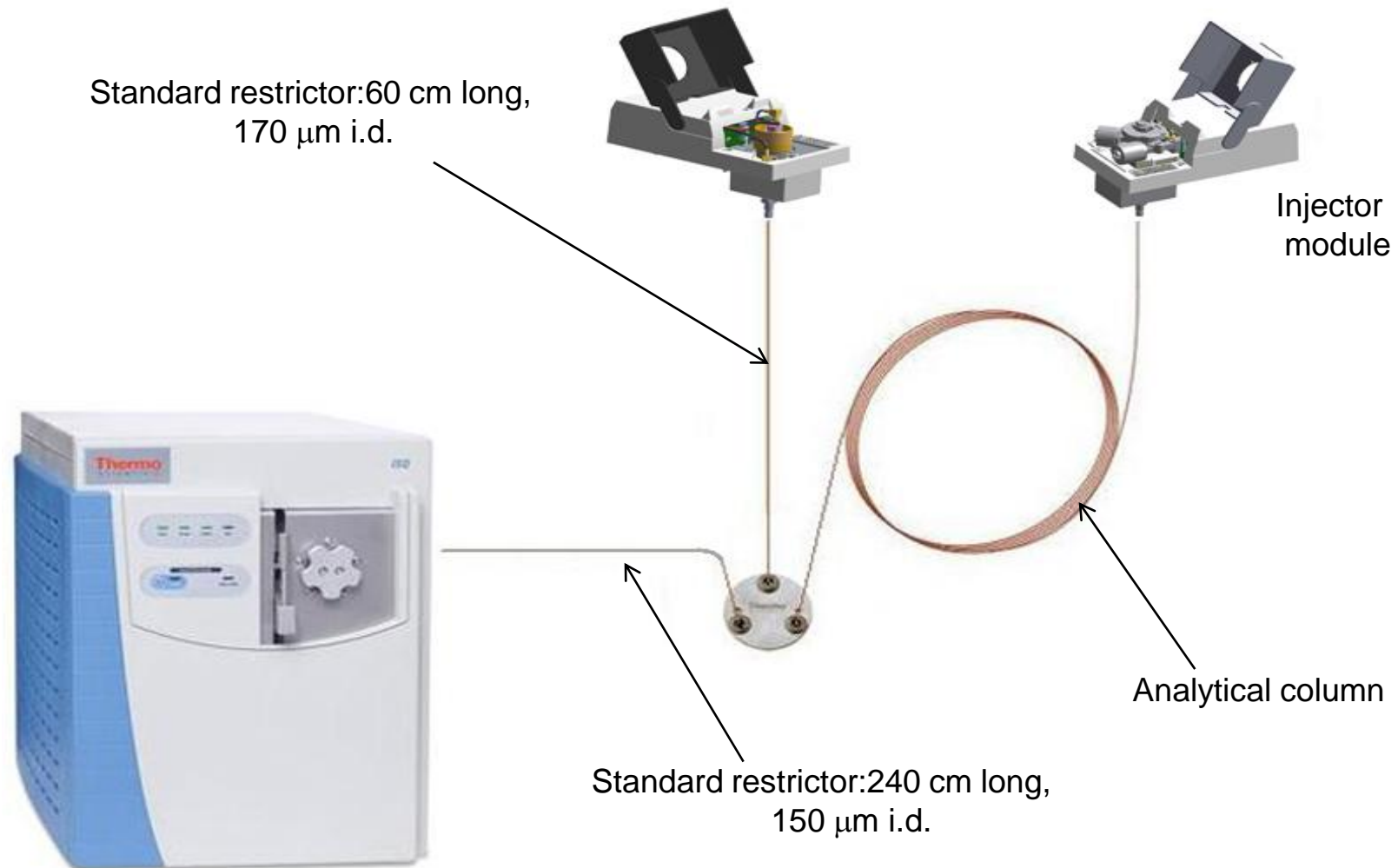



*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C

**2016 GCMS Single Quadrupole**

The world leader in serving science

# Mass Spec Components Same for All Product Offerings

Designed With Intention **ISQ Series - Intelligent Single Quad**



**ISQ QD**  
Single Quadrupole Mass Spectrometer



**ISQ QD 300**  
Single Quadrupole Mass Spectrometer




**ISQ LT**  
Single Quadrupole Mass Spectrometer

Unique Dual-Filament




ExtractaBrite™ ion


s-shaped ion guide

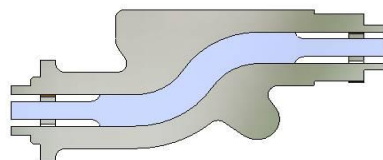

Solid, homogeneous, and maintenance-free metal quadrupole rods

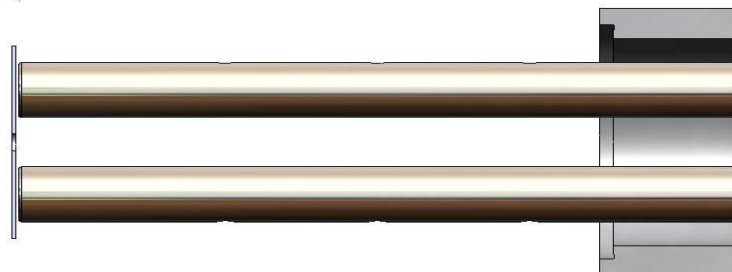

DynaMax™ XR detection

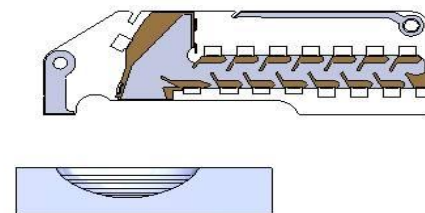

# ISQ Series Technical Specifications

Specification	ISQ QD	ISQ QD 300	ISQ LT
In <b>EI</b> mode, 1 µL of 1 pg/µL OFN m/z 272 , scanning from 50 – 300 Da (*)	<b>750:1 (*)</b>	<b>2000:1 (*)</b>	<b>2000:1 (*)</b>
In <b>PCI</b> mode, 1 µL of 100 pg/µL BZB m/z 183, scanning from 80 – 230 Da	Not Available	Not Available	300:1
In <b>NCI</b> mode, 2 µL of 100 fg/µL of OFN m/z 272, scanning from 50 – 300 Da	Not Available	Not Available	<b>2000:1</b>
IDL for m/z 272 in SIM mode in EI	<b>≤ 10 fg</b>	<b>≤ 5 fg</b>	<b>≤ 1,5 fg</b>
Ionization Modes	EI only	EI only	EI and CI
Carrier Gas Choice	He only	He or H <sub>2</sub>	He or H <sub>2</sub>
Instant Connect Helium Saver Injector	Compatible	Compatible	Compatible
Turbomolecular High Vacuum Pump	66 L/s	300 L/s	300 L/s
Vacuum Probe Interlock	Not Available (No Direct Probes)	Not Available (No Direct Probes)	VPI as standard Optional Direct Probes
Automation Features	Optional t-SIM/AutoSIM	t-SIM/AutoSIM	t-SIM/AutoSIM
Quadrupole Scanning Speed	20,000 Da/s	20,000 Da/s	20,000 Da/s

**(\*) demonstrated at installation**

# ExtractaBrite Ion Source



**Repeller**

**Patented RF Lens**

- Patented RF lens protects ion guide and quads
- RF lens is removed with source cartridge without venting
- RF lens kept hot to keep clean
- Repeller designed to overcome ion burn once it forms





# ExtractaBrite Ion Source

## ***ExtractaBrite Ion Source***

Lenses 1, 2 and 3 plus the RF lens, along with the repeller and ion volume are contained in a sleeve, all of which can be removed from **ISQ LT** while the system is under vacuum.

The ability to remove the entire ion source from **ISQ LT** without venting means you are never down and always productive.

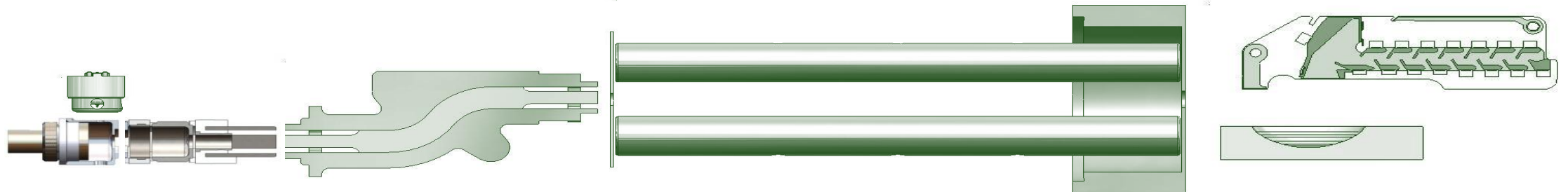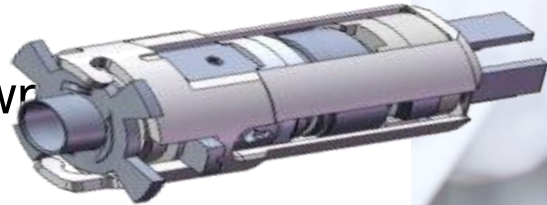

# Removable Source Under Vacuum



# GCMS Analyzer Kits



- Aim to target routine Environmental Laboratories with proven robustness, sensitivity and high sample throughput
- Compatible with ISQ QD and ISQ LT
- Ready to use package with method, troubleshooting and reporting
- Dedicated consumables and spares to complete the package
- Compatible with both Chromeleon and TraceFinder
- Compatible with Helium Saver



*Target Launch Q3*

Dedicated to **semivolatile contaminants in environmental samples**  
(EPA 8270D)



Dedicated to **semivolatile contaminants in drinking water**  
(EPA 525.3)

*Target Launch Q4*



## One Minute GC & GC-MS Highlight

### New Product: U.S. EPA 8270D Analyzer Kit

#### Update your lab and overcome the challenges of semivolatile (SVOC) analysis

Get a complete solution designed to provide seamless migration to Thermo Scientific GC-MS instrumentation for semivolatile organic compound (SVOC) analysis. The Thermo Scientific™ U.S. EPA 8270D Analyzer Kit removes the need for new method development offering a validated methodology that includes instrument and processing methods, as well as an EPA environmental specific report package. All components of the kit were developed to fulfill EPA methods requirements.

Click the links below to learn more about the new EPA 8270D Analyzer Kit:

[Link to One Minute Highlight](#)



## EPA 8270D Analyzer Kit



**Thermo Scientific™** Related applications: [Industrial Mass Spectrometry](#)

Get a complete solution designed to provide seamless migration to Thermo Scientific™ GC-MS instrumentation for semi-volatile organic compound (SVOC) analysis. The Thermo Scientific™ EPA 8270D analyzer kit removes the need for new method development offering a validated methodology that includes instrument and processing methods, as well an EPA environmental specific report package. All components of the kit were developed to fulfill EPA methods requirements.

[Contact us for support](#)

Catalog number	Product Size	Description	Software	For Use With	Price (GBP)
1R120400-8270	-	88 Validated Compounds Linearity Range: 0.2-200ppm (with same column and liner)	Chromeleon 7.2 SR4 MUB or newer; TraceFinder EFS 4.1 or newer	Thermo Scientific™ ISQ™ QD Single Quadrupole GC-MS System and Thermo Scientific™ ISQ™ LT Single Quadrupole GC-MS System	<a href="#">Request A Quote</a>

[EPA 8270D Web page](#)

EPA 8270D kit is a new, ready-to-use application solution delivering maximum up-time and more reliable results to Environmental customers



## Kit concept

- Ready to use with Method 8270D, includes:
  - Set of consumables (column, liner, septa, etc.)
  - Pre-loaded processing methods, templates and e-workflows
  - Videos with training on most common sw activities
  - App note AN10522, troubleshooting guide

## Product offering

- Available with ISQ QD and ISQ LT
- Chromeleon (SR4 Mub HF) or TraceFinder EFS compatible
- Compatible with Helium Saver
- Splitless method 0.2-50ppm and split method 2-200ppm (same column and liner)

# EPA 8270D kit p/n 1R120400-8270 BOM (Bill Of Materials)

Part number	Description	GRP, \$USD
26098-1425	Column, TG-5MS w/ 5m guard,30m X 0.25mm X 0.25µm	\$637,00
453A1925	Injection liner,4mm ID x 105 mm, pack of 5	\$109,00
290VA191	Graphite/Vespel Ferrule, 0.1 - 0.25 mm, pack of 10	\$48,00
31303233	Septa BTO Coated, 11 mm, pack of 50	\$81,00
1R120400-SCA	System Conditioning Agent	\$85,41
1R120479-0060	EPA 8270D CD (contains methods, user guide, video training, app note, troubleshooting guide) compatible with TF and CM	\$500,00
	Total Reference Price <u>of all parts</u>	\$1460.41
<b>1R120400-8270</b>	<b>EPA 8270D Kit Global Reference Price</b>	<b>\$800,00</b>



- TraceFinder methods and templates
- Chromeleon e-Workflow (with methods and reports)
- EPA 8270D user guide
- EPA 8270D application note
- Tutorial Video



thermo scientific

APPLICATION NOTE

## Optimizing the Analysis of Semi-volatiles by EPA Method 8270D

No. 10522

**Authors:** David Steiniger, Tommaso Albertini, Richard Law

**Keywords:** GC-MS, Semi-volatiles, Helium Saver, ISQ, 8270

### Introduction

The U.S. Environmental Protection Agency (EPA) released the first Semi-Volatile Organic Compounds method by Gas Chromatography/Mass Spectrometry (EPA 8270) at the end of 1980 and is a common method used in almost all environmental labs. Since then, single quadrupole mass spectrometers have become thousands of times more sensitive and the relative abundances have changed.

### Experimental Conditions

Thermo Scientific ISQ Single-Quadrupole GC-MS System


The method was tested on 5 different ISQs units to assess method transferability and instrument-to-instrument variability. Both ranges (0.2–50 ppm and 2–200 ppm) validated using Thermo Scientific Instant Connect Helium Saver Module (P/N 19070013) and Instant Connect Helium Split-Splitless SSL module (P/N 19070010). The liner is Injection line, 4 mm ID x 105 mm (P/N 453A1925). ISQ was operated in full-scan mode and the Thermo



## EPA 8270D Kit for ISQ Series Mass Spectrometers

User Guide

1R120555-0010 Revision A July 2016

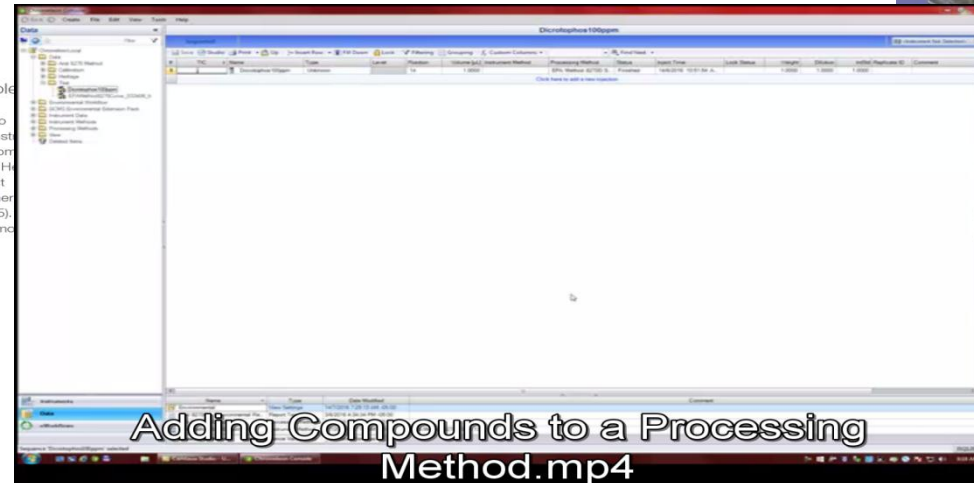

thermo scientific



# How to order EPA 8270D kit

1°

- Select GC mainframe for GC/MS (Trace 1300 or Trace 1310)

2°

- Optional stage, select auto sampler

3°

- Select GC injector module (SSL or He Saver)

4°

- Select SW (TraceFinder EFS or Chromeleon) and Computer

5°

- Select MS detector (ISQ QD or QD 300 or LT)

6°

- Select EPA 8270D kit (p/n **1R120400-8270**)

7°

- If needed or requested select application support at customer site

[LINK TO SALES TRAINING WEBINAR](#)



- Integrated control for TSQ8000Evo
  - Performance Improvements with Large Datasets
- Pesticides Analyzer
  - Update to GCBU Pesticides Analyzer DVD
  - Will include method templates, report templates, compound database and workflow instructions
- GCMS EPA Water Analysis
  - Environmental Extension Pack containing eWorkFlows, method templates and report templates, view settings, custom variables
- Deconvolution SW for unknowns and RI
  - Target release in Q4 2016




*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C

2016 GCMSMS TSQ 8000 EVO , TSQ 8000 DUO

# Triple Quadrupole MS – Principle of Operation

## Quantitation of target compounds in matrix samples

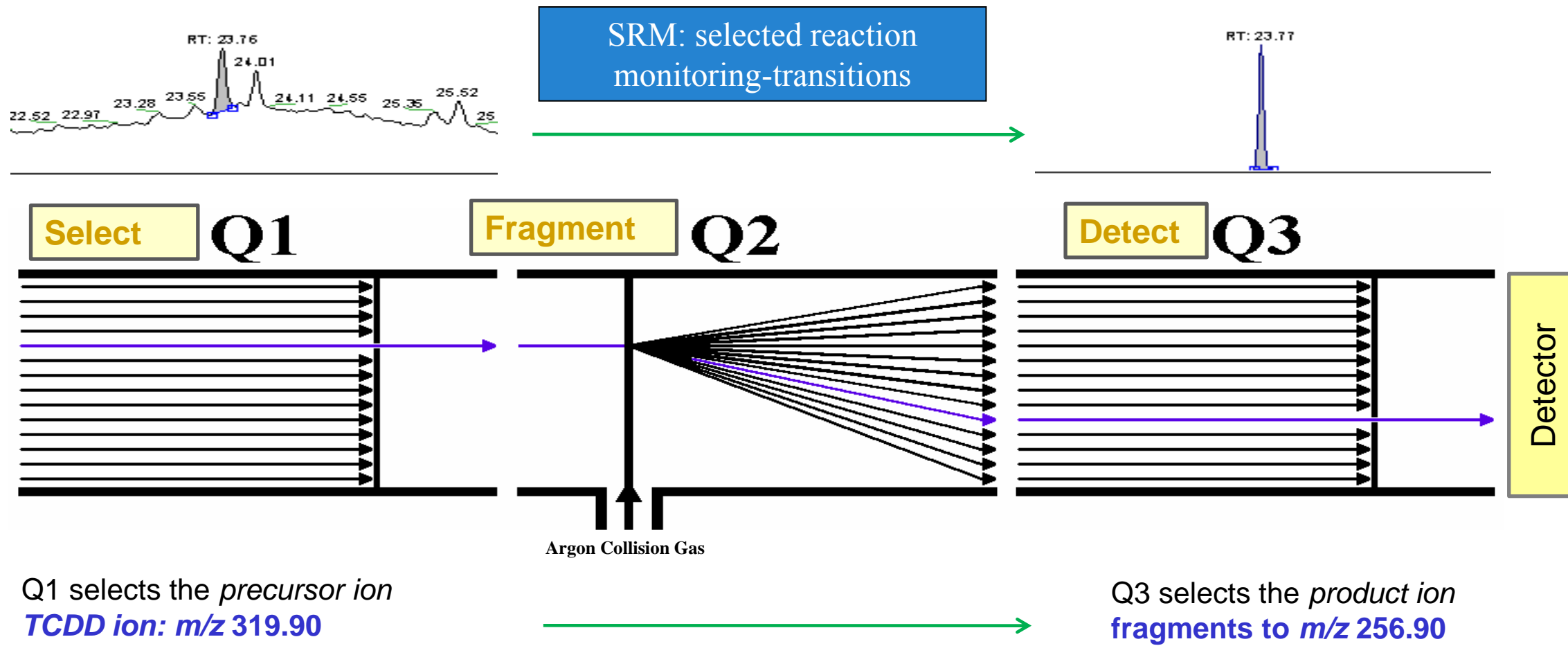

# GCMSMS Analyzer Kits



- Aim to target routine Environmental and Food Safety Laboratories with proven robustness, sensitivity and high sample throughput
- Compatible with TSQ DUO and TSQ8000EVO
- Ready to use package with method templates, report templates, compound database and workflow instructions
- Dedicated consumables and spares to complete the package
- Compatible with both Chromeleon and TraceFinder
- Compatible with Helium Saver



## Pesticide Analyzer

CDB extended to 1001 compounds with 3 optimized transitions


## Dioxins Analyzer

Dedicated workflow and SOP for Dioxins Analysis




*Target Launch Q3*

# Drugs of Abuse Testing Laboratory Workflow

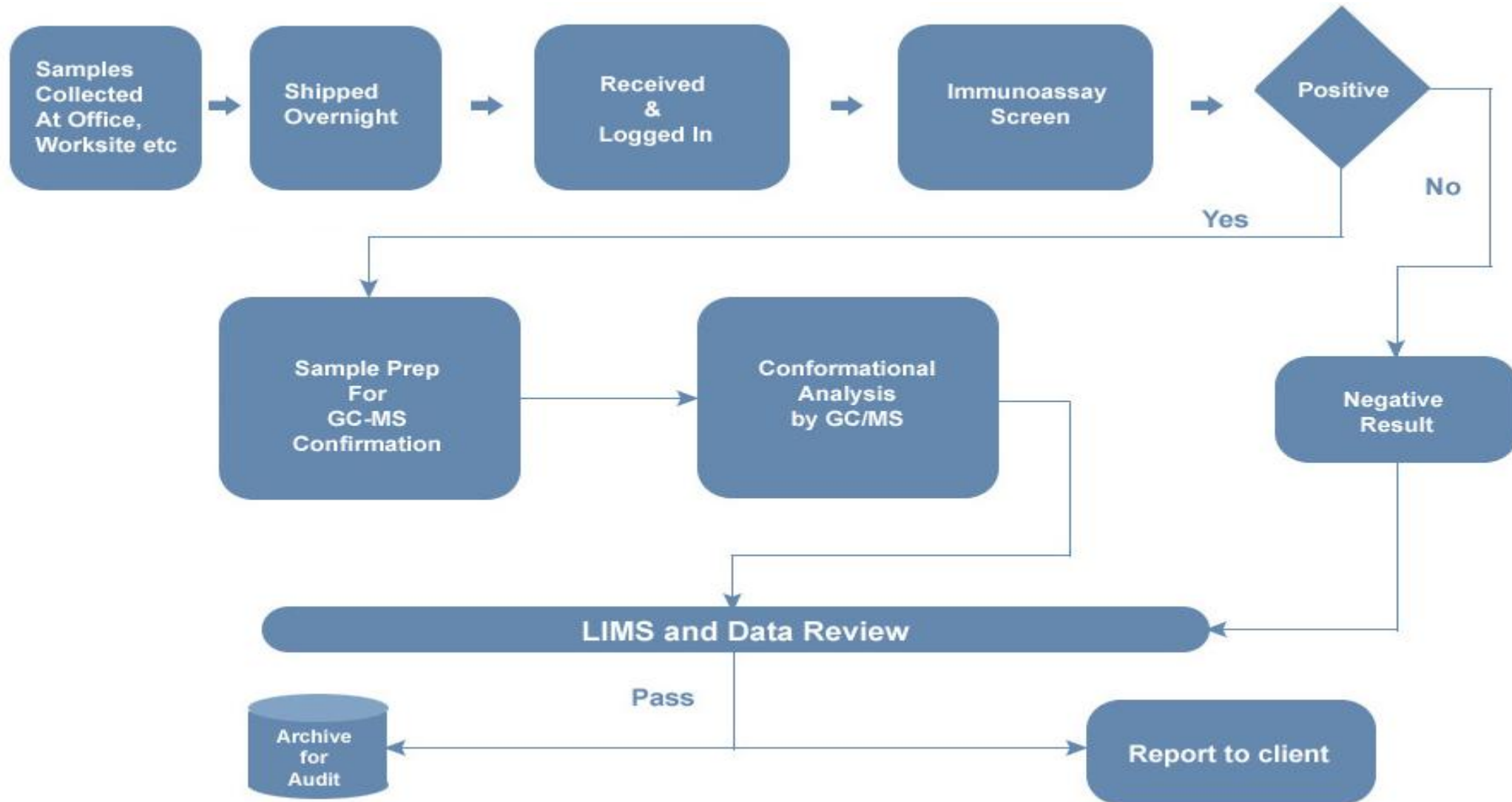

# Sample prep: injecting blood, sweat and tears on a GC

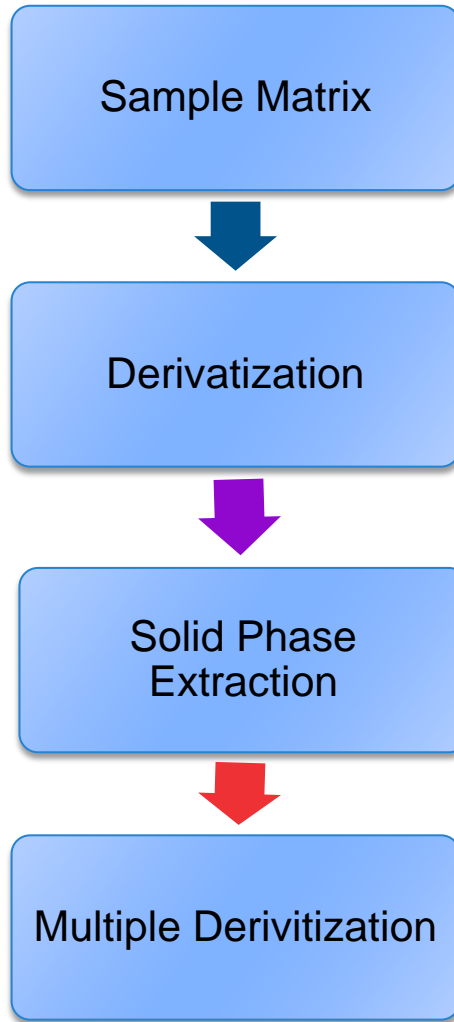

- Sample preparation reagents and solutions for various matrices
  - Blood, Urine, Oral Fluid hair and other bodily fluids
- Accurate separation of components analogues
  - Avoid misinterpretation of target compounds
  - Avoid metabolic product exclusion or interference
    - Opioids - Keto-moeity derivitization
- Removal of noise associated with sample matrix components
  - Higher recovery, Better reproducibility
  - Less noise = Better chromatography
    - Hyper Sep Verify - CX
- Enhance target compound stability
  - Better detection

- **Why test THCA?**

Delta 9 THC (THC) is the most commonly abused illicit drug.

The metabolite THCA indicates the ingestion of THC as opposed to passive contamination of the hair through the air.

- **Why test Hair?**

Long detection time windows

*Urine/blood/oral fluid: minutes to days detection window*

*Hair: months to years detection window*

- **Applications for Hair testing**

Death investigations

Workplace drug testing

Drug-facilitated crimes


Violation of probation or parole

Drivers license re-instatement




# Determination THCA in Hair

- Delta 9 THC (THC) is the most commonly abused illicit drug.
- The metabolite THCA indicates the ingestion of THC as opposed to passive contamination of the hair through the air.
- Hair is a very complex and variable matrix even after sample cleanup

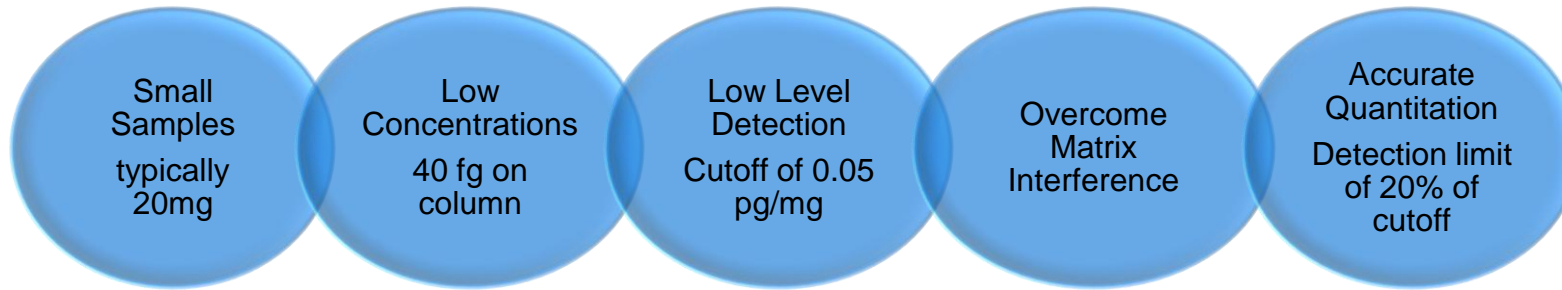

<b>Agency/Organization</b>	<b>Cutoff Level</b>
Society of Hair Testing (SoHT)	0.20 pg/mg
Gesellschaft für Forensische und Toxikologische Chemie (GFTCh)	0.05 pg/mg
Substance Abuse and Mental Health Services Administration (SAMHSA)	0.05 pg/mg

# Determination THCA in Hair

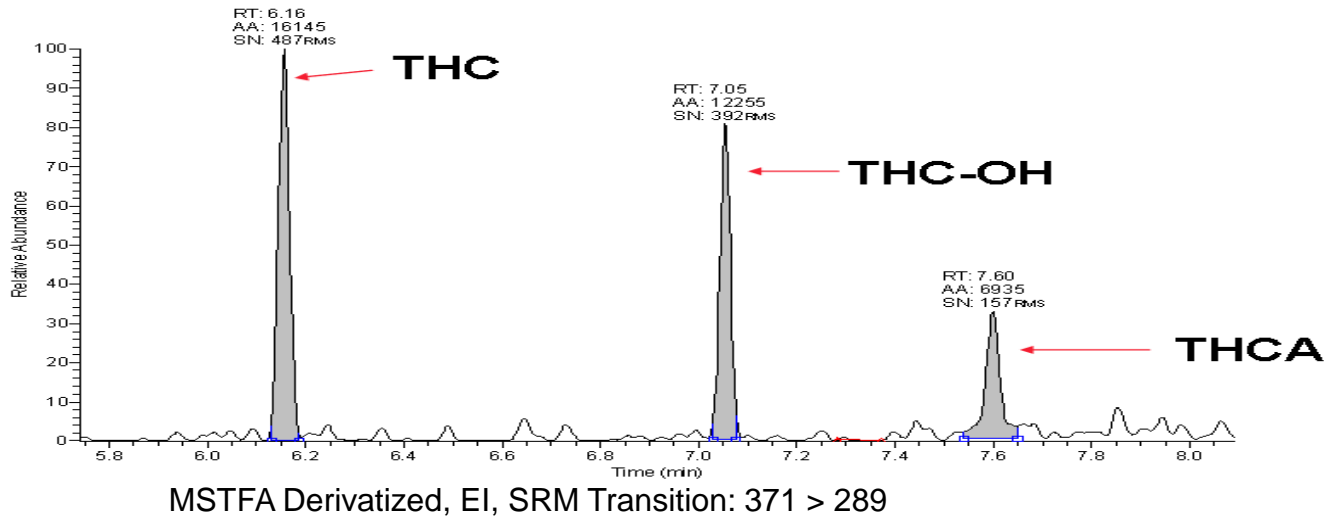


- Sensitivity of GC TSQ 8000 EVO to THCA
- On-Column Injection of 8fg
- Use of a 0.05 pg/mg cutoff
- Demonstrates the sensitivity sufficient for published cutoffs

- **NCI** - Specific transitions can be defined in NCI (Methane used as reagent gas), enhancing selectivity for both THC-COOH and the corresponding ISTD D3

Compound	Precursor ion (m/z)	Daughter ion (m/z)	Coll cell (mtorr)	Width	Time	CE	Q1 PW	Q3 PW
THC-COOH	620.2	493	1.5	0.5	0.05	15	0.7	0.7
THC-COOH	620.2	532	1.5	0.5	0.05	15	0.7	0.7
THC-COOH d3	623.2	496	1.5	0.5	0.05	15	0.7	0.7
THC-COOH d3	623.2	535	1.5	0.5	0.05	15	0.7	0.7

# THCA in Hair - NCI transitions

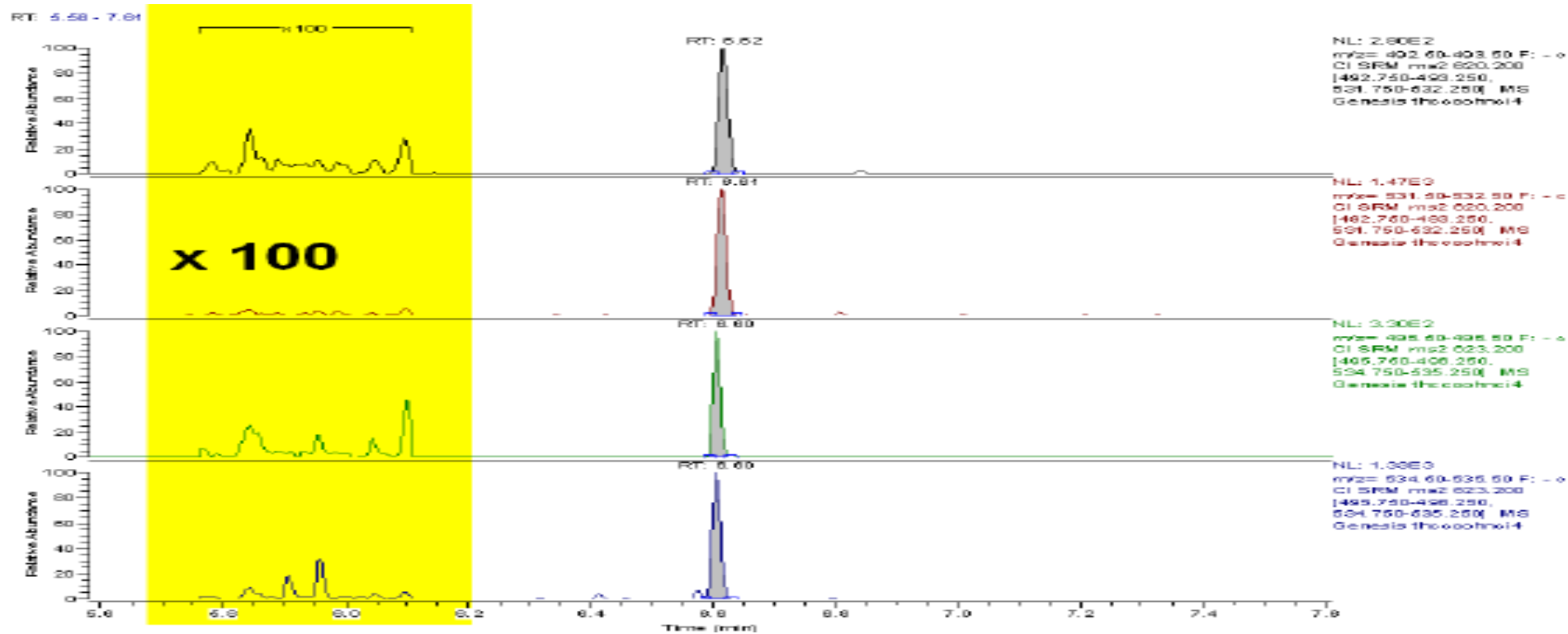

PFAA Derivatized, Negative Ion Chemical Ionization - **50 fg** THCA and ISTD-D3

- Demonstrates very low noise level relative to the signal
- Zoom in the baseline by 100 times before we start seeing significant noise
- Possible to drop the concentration by a factor of 10 and still achieve the selectivity to accurately quantitate THCA



*Unleash the*  
**value!**

**ThermoFisher**  
S C I E N T I F I C

**2016 Qexactive GC & Exactive GC Orbitraps**

# Exactive GC / Q Exactive GC



Redefining Routine GC-MS

Thermo Scientific Exactive GC  
Orbitrap GC-MS System

**50,000 maximum Resolution**

**18 Hz maximum speed**

Mass Accuracy

<1ppm RMS for Internal Calibration

< 3ppm RMS for External Calibration




**100,000 maximum Resolution**

**18 Hz maximum speed**

Mass Accuracy

<1ppm RMS for Internal Calibration

< 3ppm RMS for External Calibration

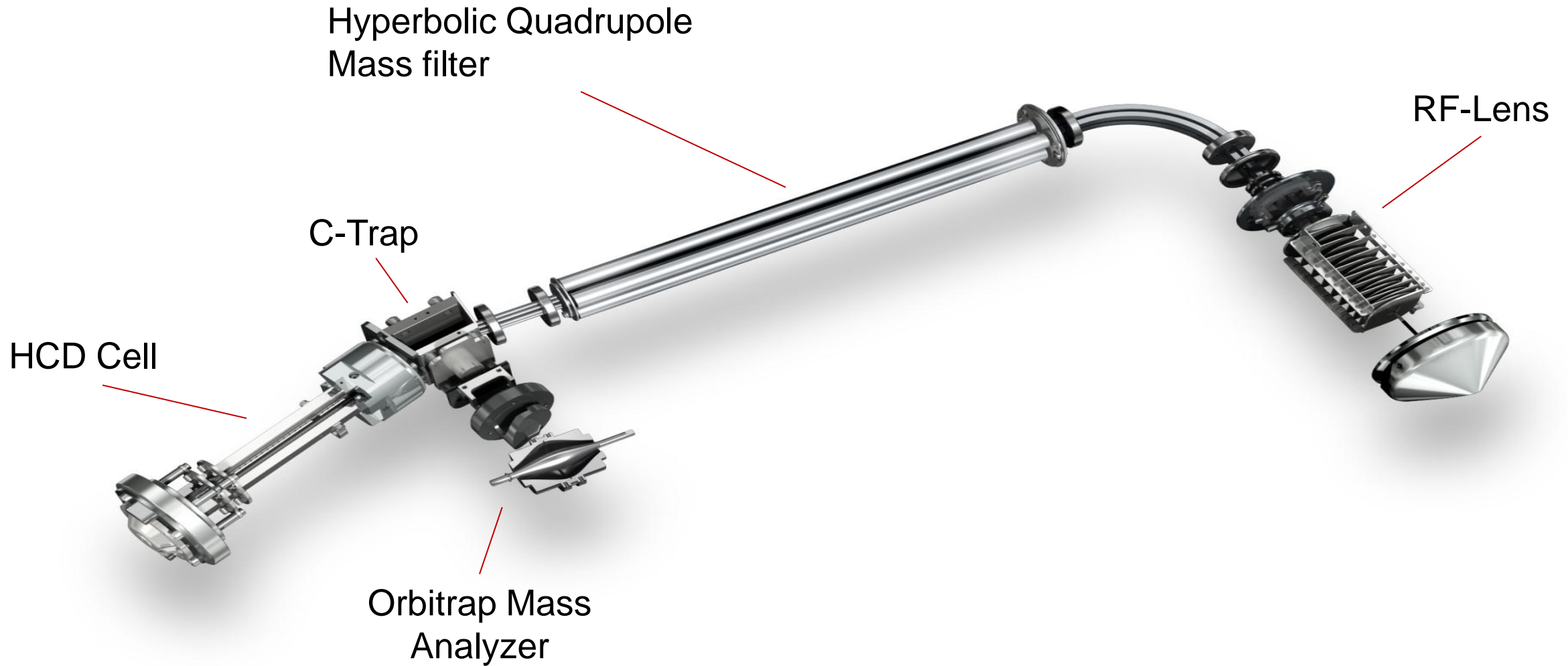


# Exactive GC Internal Components




Note for reasons of internal product positioning, AQT quadrupole won't be highlighted in marketing material

# Instrument Overview

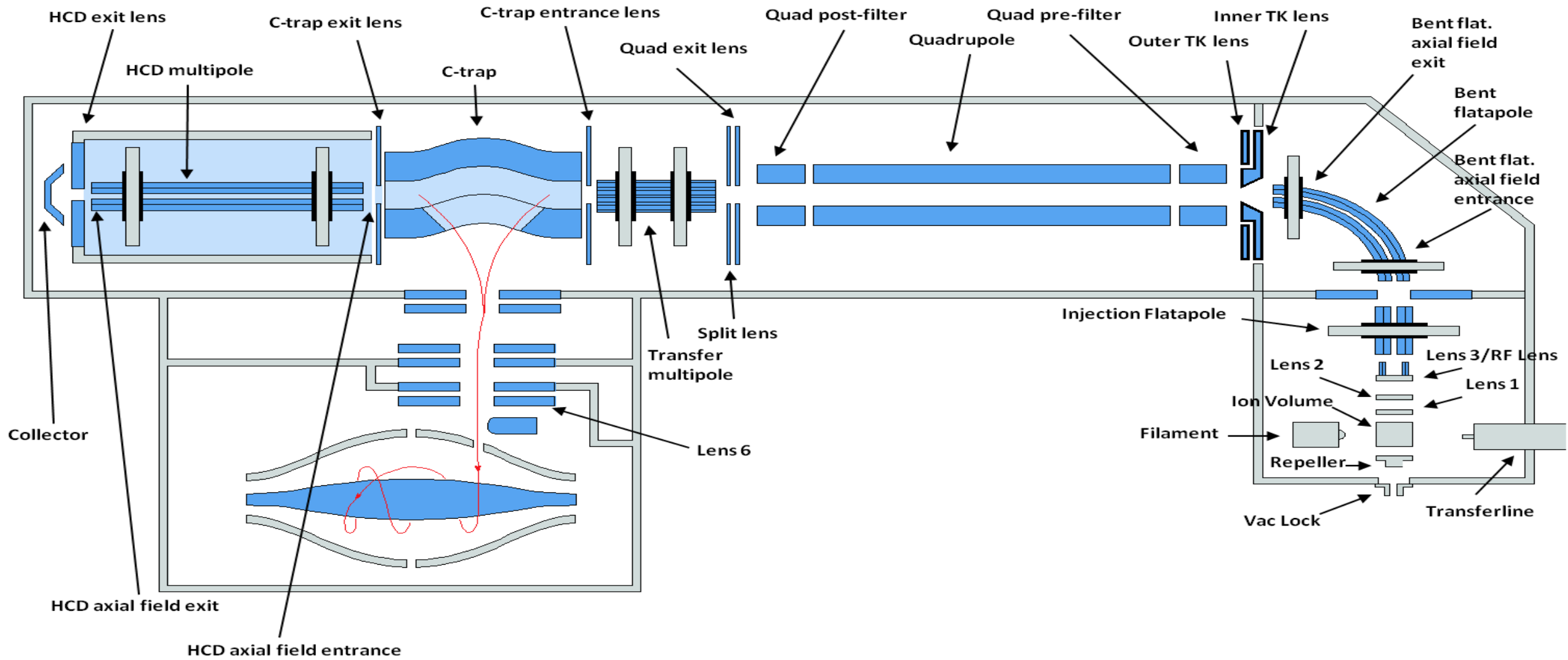




# 90° Bent Flatapole

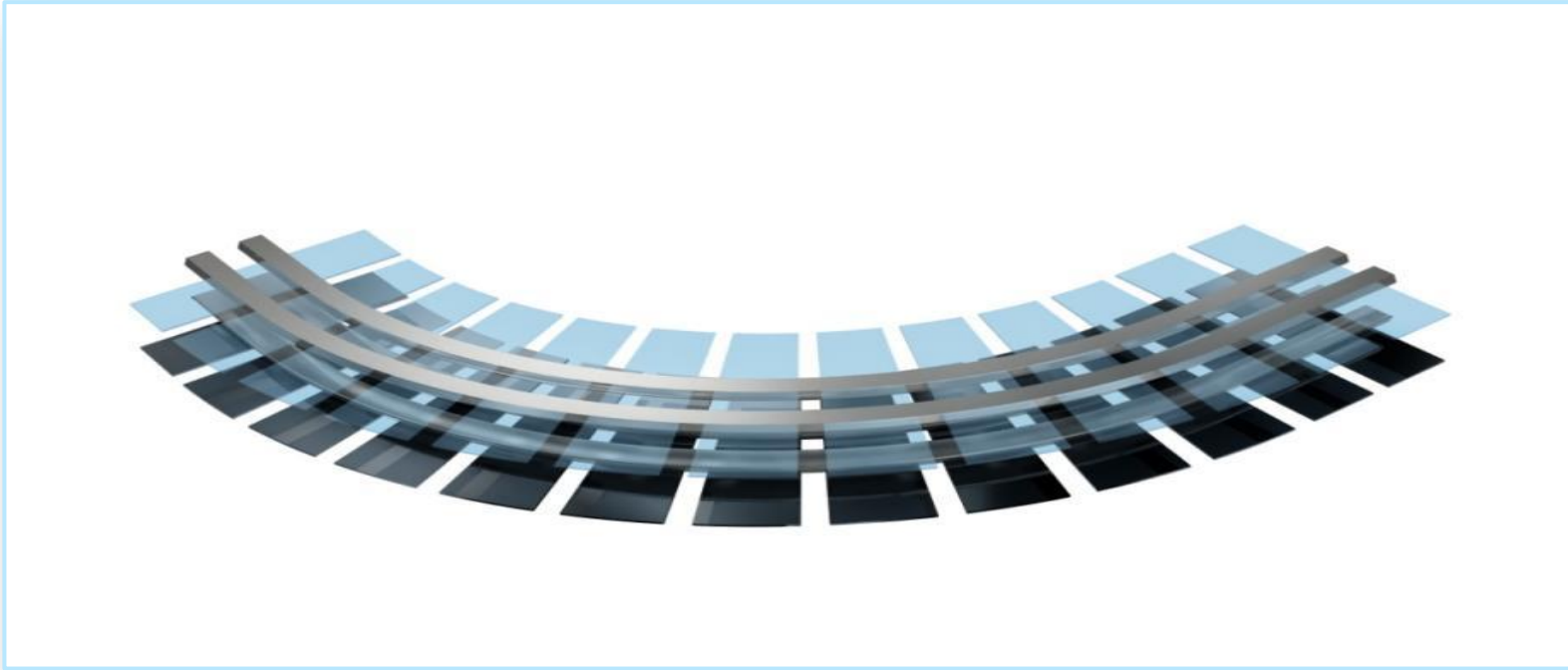


- Similar design to Q Exactive Plus and Q Exactive HF
- Axial field to optimize ion transfer
- Neutrals are removed from ion beam

# Video of Basic Orbitrap Scan Modes

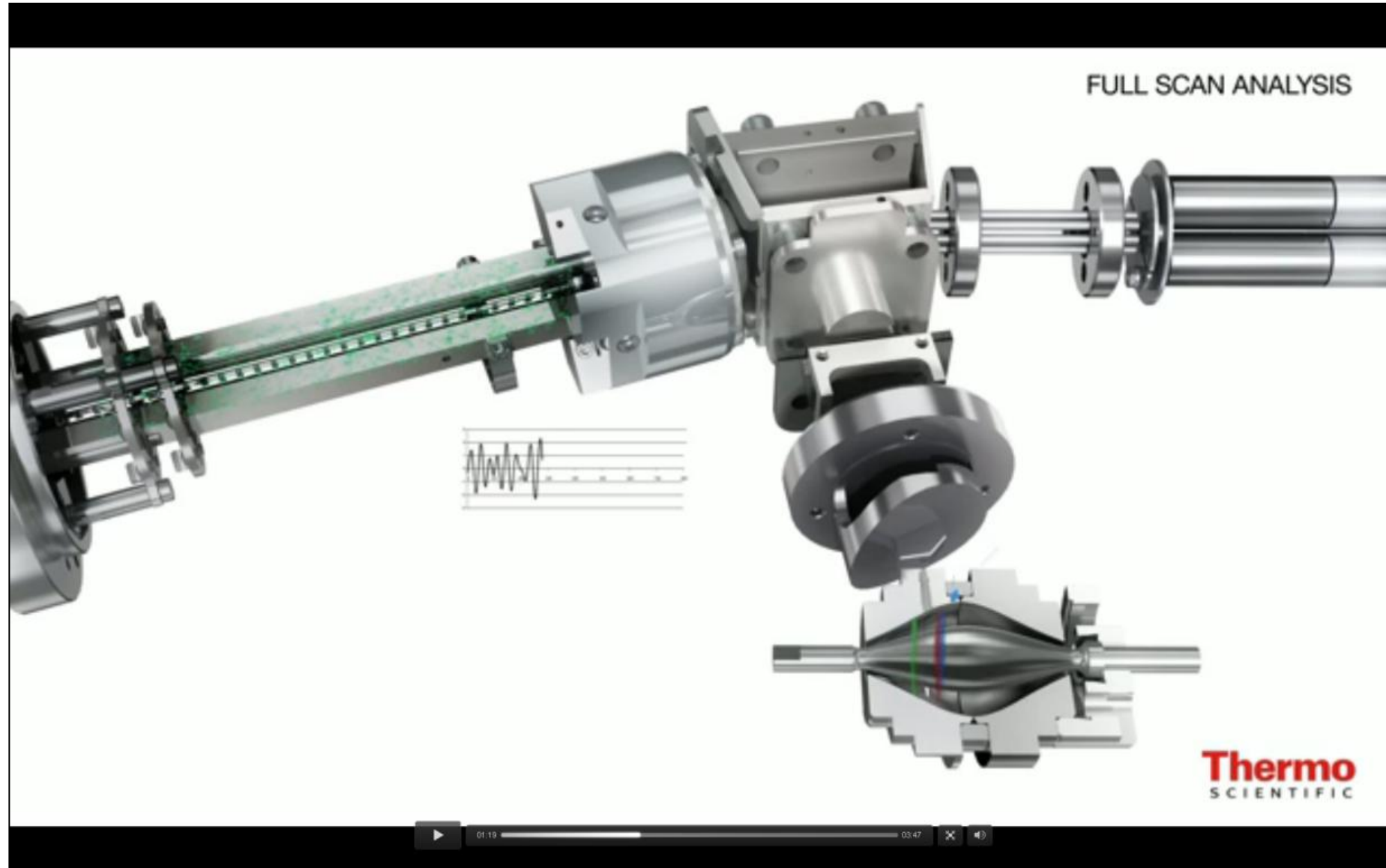

# Retaining the Ions in the Orbitrap

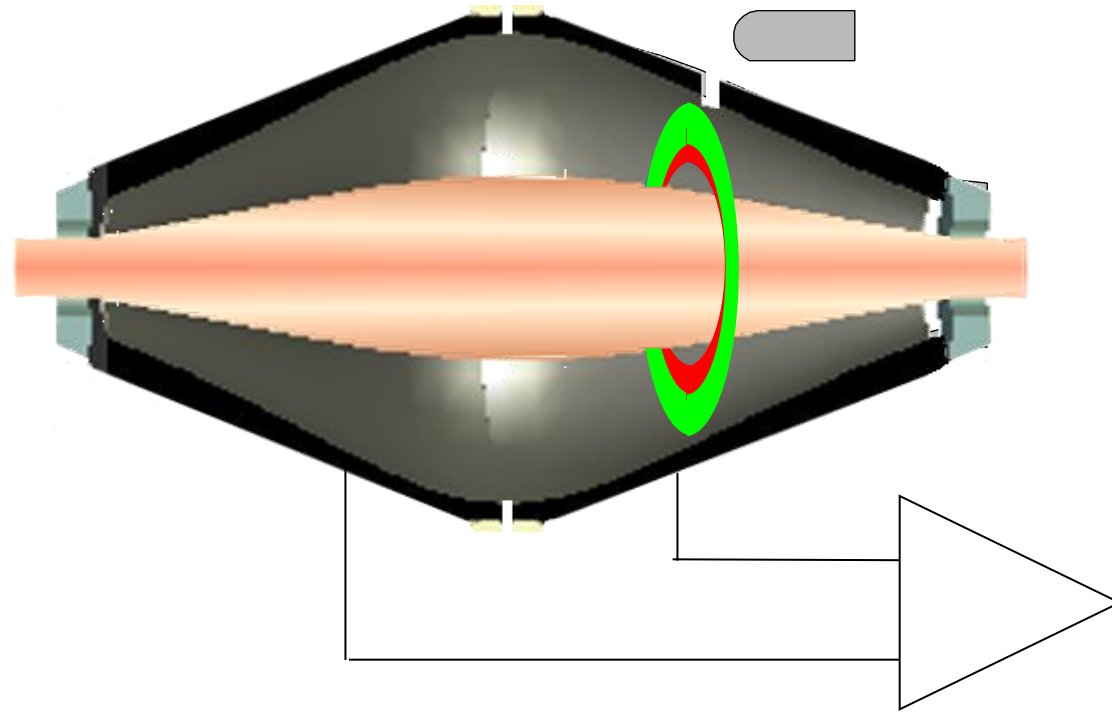


$$\omega = \sqrt{\frac{k}{m/z}}$$

- Frequency of axial oscillations are independent of initial conditions of ions entering trap
- Therefore these oscillations used for mass determination
- Image current measured by outer split electrodes (no electron multiplier to replace!)
- Ion frequencies determine by complex superposition of measured ring oscillations through Fourier transform

# What the Orbitrap Looks like

Inner-electrode

Split outer-electrode

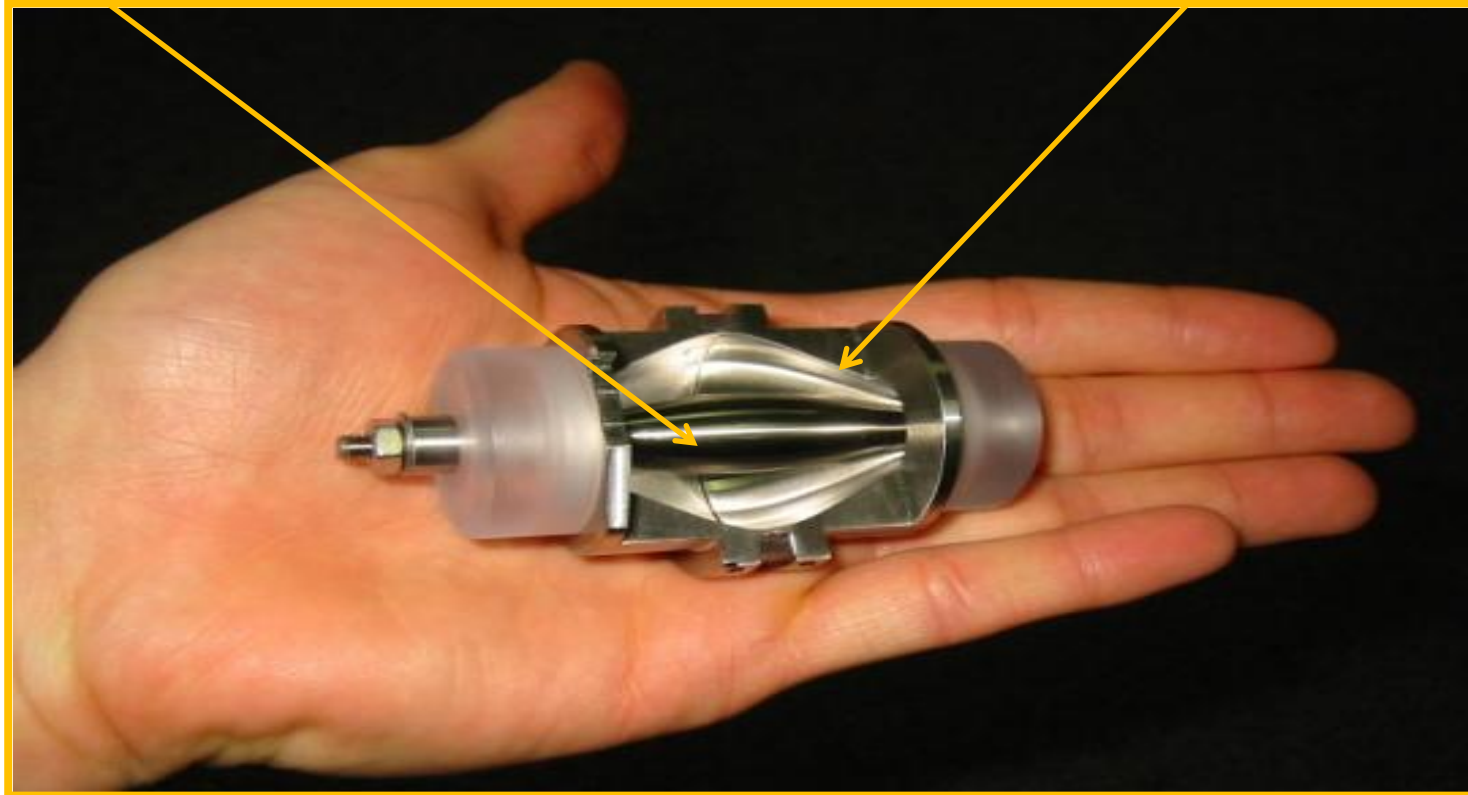

# Definitions to remember

- **Nominal mass:** the mass of an ion calculated using the *integer mass of the most abundant isotope of each element*.
  - ignores the mass defect, where H=1, C=12, O=16
- **Monoisotopic mass:** the mass of an ion calculated using the *exact mass of the most abundant isotope of each element*.
  - includes the mass defect, where  $^1\text{H}=1.0078$ ,  $^{12}\text{C}=12.0000$ ,  $^{16}\text{O}=15.9949$
- **Average mass:** the mass of an ion calculated using the *relative average isotopic mass of each element*.
  - where  $\text{C}=12.0111(12.000000) (0.9890) + (13.003355) (0.0110) = 12.011 \text{ Da}$
  - H=1.00797, O=15.9994
- **Isotopic abundance:** the naturally occurring distribution of the same element with different atomic mass.
  - ex:  $^{12}\text{C}=12.0000 = 98.9\%$   
 $^{13}\text{C}=13.0034 = 1.1\%$

# How accurate is your mass?

- Mass accuracy

$$\Delta m / z = \frac{m_{meas} - m_{true}}{m_{true}} \cdot 10^6$$

- Quadrupole MS

$$\Delta m / z = \frac{500.1 - 500.0}{500} \cdot 10^6 = 200 \text{ ppm}$$

- Orbitrap MS




$$\Delta m / z = \frac{500.10314 - 500.10214}{500.10314} \cdot 10^6 = 2 \text{ ppm}$$

- Typical mass accuracy capability for various MS types:

<b>Type</b>	<b>Mass Accuracy</b>
FT-ICR-MS	0.1 - 1 ppm
Orbitrap	0.5 - 1 ppm
Magnetic Sector	1 - 2 ppm
TOF-MS	3 - 5 ppm
Q-TOF	3 - 5 ppm




*Source: Metabolomics Fiehn's lab*

# Compound ID confirmation

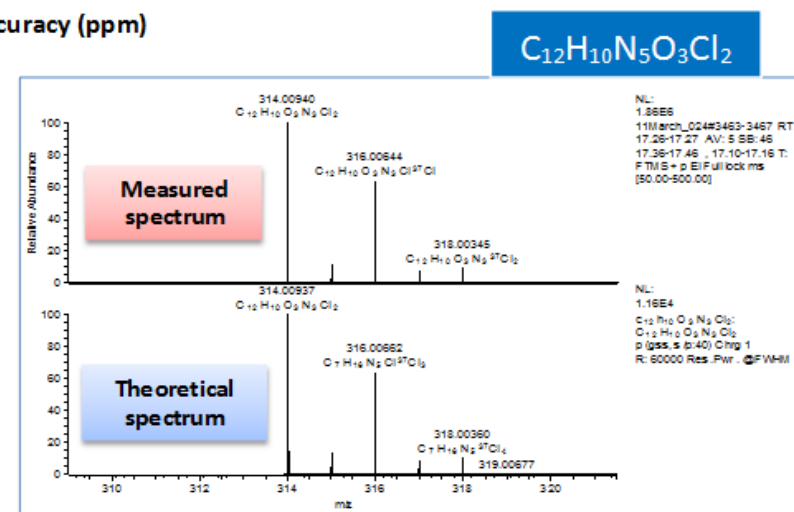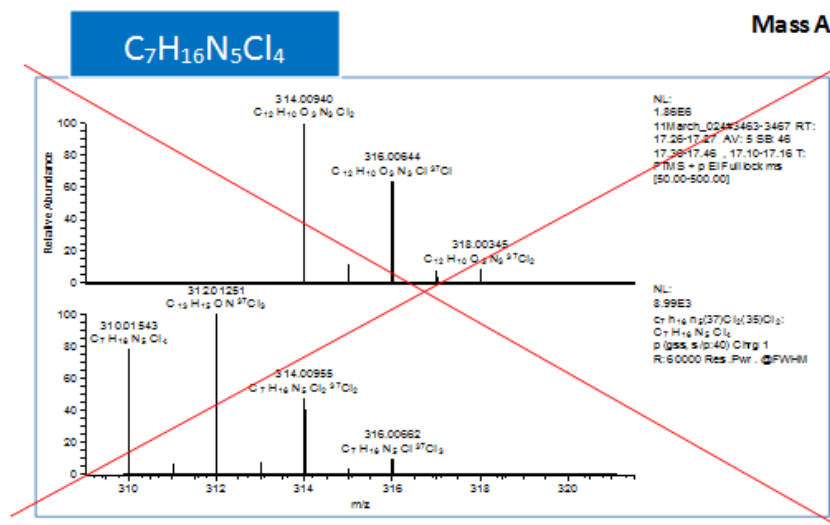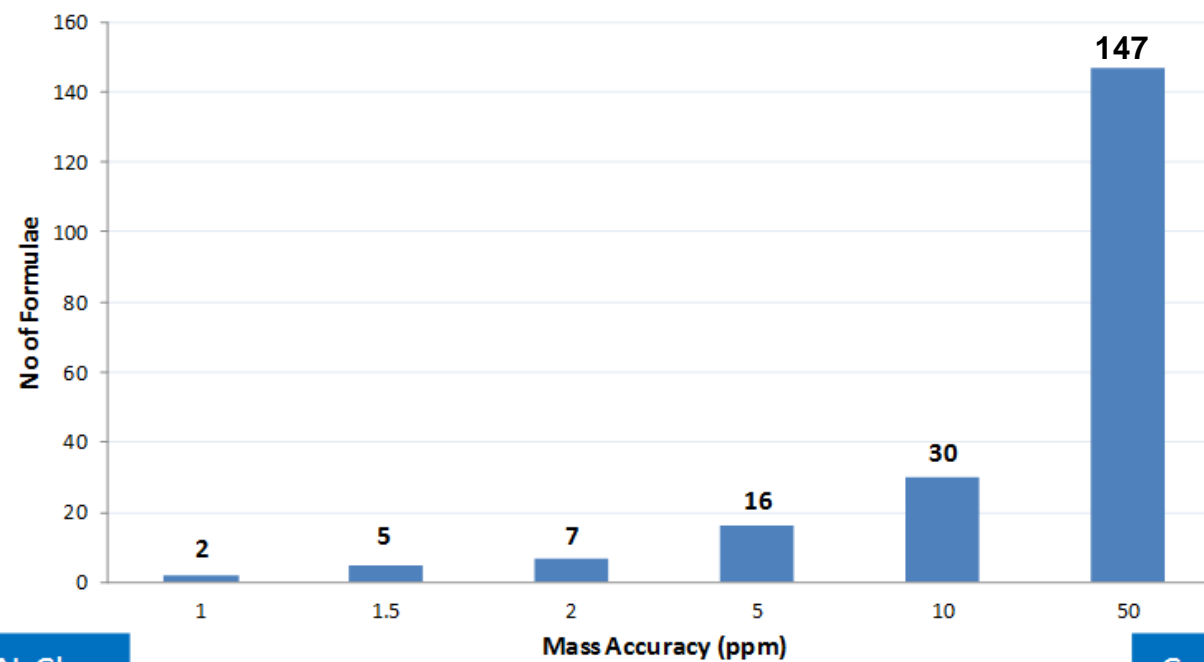



# The Need for High Resolution

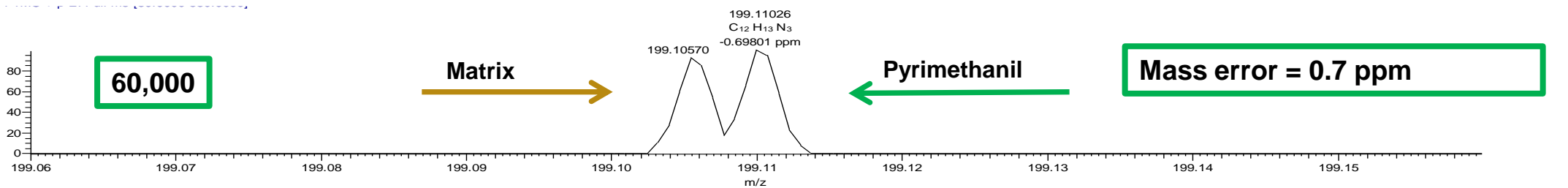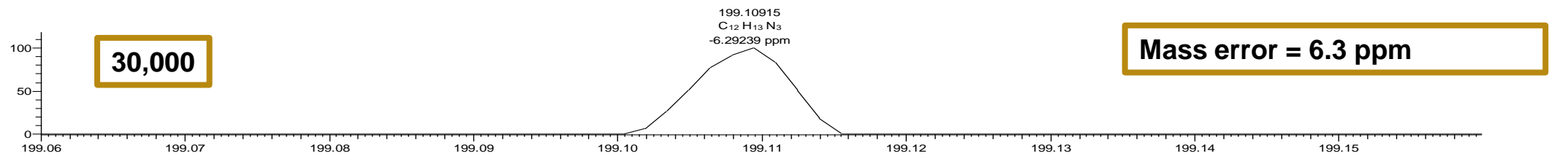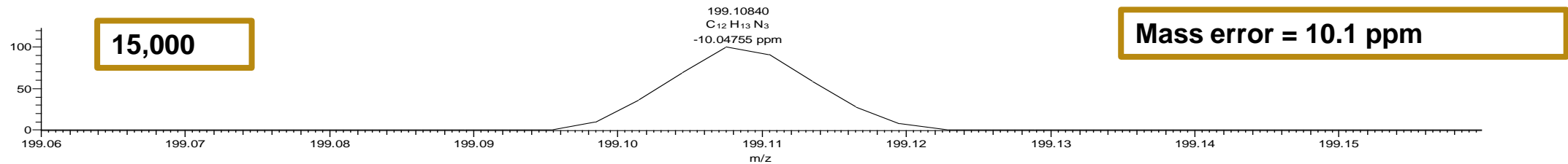


## Thermo Scientific Q Exactive GC Hybrid Quadrupole-Orbitrap GC-MS/MS System

Unprecedented Depth in GC-MS Analysis



- Thank you for your attention!
- Q&A

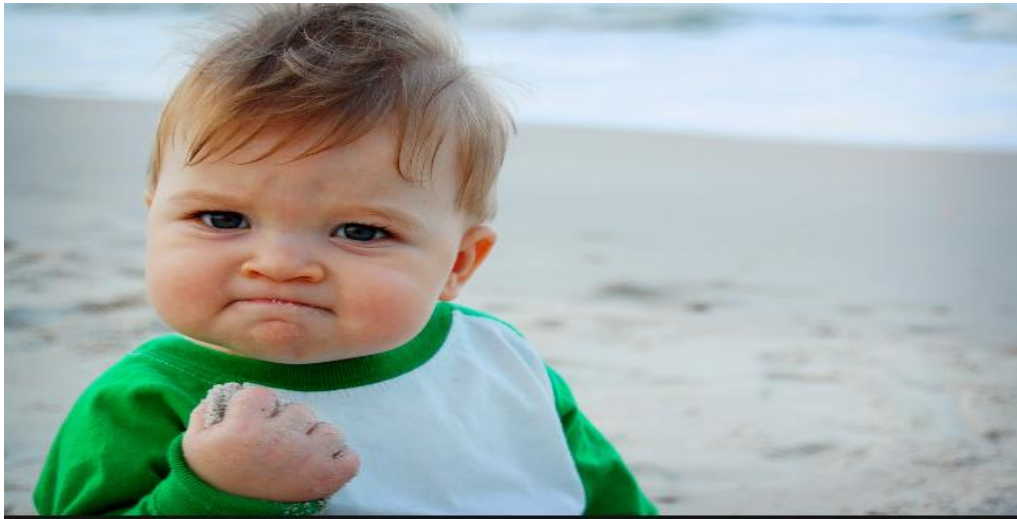

***We are changing the game!***

