

PURELAB flex

Pure and ultrapure laboratory water from a single system





Your laboratory water specialists

ELGA has been a trusted name in water purification for over 50 years, pioneering innovative technologies and award winning product design for our customers. Part of the world's leading water services company Veolia Water Solutions and Technologies, ELGA have the most comprehensive range of water purification systems available today for laboratory research, healthcare and diagnostics.

Water is the most common reagent used in today's laboratories for research and testing applications and we understand how important accurate consistent test results are. The PURELAB flex range is designed around your needs, delivering accuracy, flexibility and ease of use with an innovative and ergonomic design.

Find a PURELAB flex to match your application



PURELAB flex 2

PURELAB flex 3&4

PURELAB flex 1





PURELAB flex 3&4

PURELAB flex 1

Daily volume	>10 Liters	<10 Liters	>10 Liters	<10 Liters	>10 Liters
Water quality	Type I ultrapure water		Type II/Type II⁺ pure water	Type III RO water	
Feedwater source	Pre-purified water	Pre-purified or tap water	Pre-purified water	Pre-purified or tap water	Pre-purified water
Typical applications	Electrochemistry Electrophoresis GFAAS (Graphite Fur Absorption Spectrop HPLC IC (Ion Chromatogra ICPAES (Inductively C Atomic Emission Spe ICPMS (Inductively C Spectrometry) Mammalian and bac Molecular biology Plant tissue culture Qualitative analyses	hotometry) Joupled Plasma (ctrometry) oupled Plasma Mass terial cell culture	AAS (Atomic Absorption Spectrophotometry) Buffer and media preparation Electrophysiology FAAS (Flame Atomic Absorption Spectrophotometry) Feed to ultrapure water systems Glassware washing/rinsing General chemistry Histology Microbiological analysis RIA (Radioimmunoassay) / ELISA (Enzyme Linked Immunoabsorbant Assay) Sample dilution and reagent preparation Spectrophotometry	Autoclave feed Feed to ultrapure war Hydroponics Plant growth cabinet Stability chambers Steam generators Sterilizer feed	

Water purification made easy

Intuitive to use, ergonomic handset design

Clear water purity display for absolute confidence as you dispense

Realtime TOC monitoring for critical analytical applications

> Prioritized information shown at all times (system status, TOC, alarm)

Flexible dispensing in four different ways:

- Variable flow drop by drop or up to 2 liters per minute
- Autovolume dispense from 50ml to 60 liters and repeat dispensing
- Hands free with optional foot pedal
- Locked dispense for glassware filling

POU filters for multiple applications

- Remove endotoxins, DNase, RNase and bacteria
- Conforms to international water standards e.g. CLSI, CLRW, ISO 3696: Grade 1,2,3, ASTM D1193-06, Pharmacopeia USP, EP and JP

Optional accessories

- Wall mounting bracket
- Foot pedal
- Thermal printer
- Point of Use filters
- Leak sensor

Designed for today's laboratories



Quick and easy sanitization


Minimize microbial growth to give you absolute confidence in the water purity. The simple procedure involves one single sanitization pack so no harsh chemicals need to be ordered, handled or mixed. No additional consumables need to be replaced during this short process.


Global network, local support

Our highly trained service engineers will apply their expertise to the installation, validation and maintenance of your water purification system. We aim to provide a first-time-fix service.

Inside the PURELAB flex

PURELAB flex 2 Process Flow Diagram

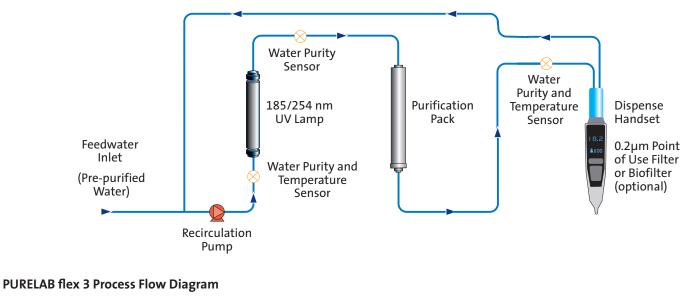

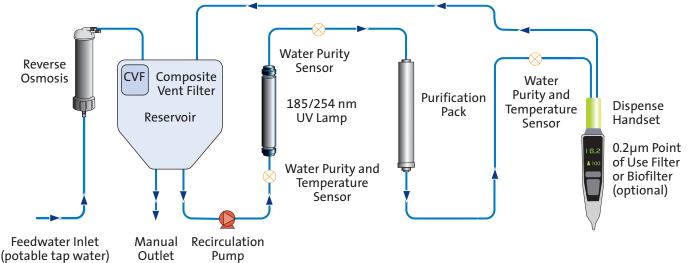

Want to know more?

Find out more or book your demonstration at www.purelabflex.com

Commitment to quality

PURELAB flex systems are supplied with a Certificate of Conformity, ensuring they have been manufactured and tested at ELGA LabWater Global Operations, whose Quality Management System is approved by Lloyds Register Quality Assurance (LRQA) and complies with BS EN ISO 9001:2008.

ELGA's manufacturing facility has been approved by LRQA and complies with the environmental management system ISO 14001:2004.

The PURELAB flex has been designed to conform to the requirements of the Waste Electrical & Electronic Equipment (WEEE) Directive and the Restriction of Hazardous Substance (RoHS) Directive.

Assured safety

The PURELAB flex is compliant with CE directives associated with safety and electromagnetic compatibility and is approved and tested by an accredited external company. The PURELAB flex is also compliant with the requirements of IEC/EN61010-1 and UL61010-1 and is marked with an ETL label showing compliance.

Treated Water Specifications

MODEL	PURELAB flex 1	PURELAB flex 1 With Purification Pack	PURELAB flex 2	PURELAB flex 3	PURELAB flex 4
Daily volume	>10 liters	>10 liters	>10 liters	<10 liters	<10 liters
Delivery flow rate – maximum	Up to 2 l/min	Up to 2 l/min	Up to 2 l/min	Up to 2 l/min	Up to 2 l/min
Inorganics (resistivity @25°C)	As per feedwater	18.2 MΩ-cm	18.2 MΩ-cm	18.2 MΩ-cm	18.2 MΩ-cm
Organics (TOC) – typical	Dependant on feedwater	Dependant on feedwater	<5ppb ¹	<5ppb ¹	<5ppb ¹
Direct from internal reservoir	N/A	N/A	N/A	Type III / RO water ²	Type III / RO water ²
Bacteria – typical (when fitted with POU Filter)	<1 CFU/10ml	<1 CFU/10ml	<1 CFU/10ml	<1 CFU/10ml	<1 CFU/10ml
Bacteria – typical (when fitted with Biofilter)	<1 CFU/10ml	<1 CFU/10ml	<1 CFU/10ml	<1 CFU/10ml	<1 CFU/10ml
Endotoxin (when fitted with Biofilter)	N/A	<0.001 EU/ml	<0.001 EU/ml	<0.001 EU/ml	<0.001 EU/ml
DNase (when fitted with Biofilter)	N/A	N/A	<20 pg/ml	<20 pg/ml	<20 pg/ml
RNase (when fitted with Biofilter)	N/A	N/A	<0.002 ng/ml	<0.002 ng/ml	<0.002 ng/ml

¹ Dependant on feedwater ² (Equal to or >) Type III RO Water

Dimensions and weights

Dimensions	Width 236mm, Depth 374mm, Height minimum 900mm, Height maximum 1020mm			Width 236mm, Depth 470mm, Height minimum 900mm, Height maximum 1020mm	
Operational weight	10kg (22lb)	10.5kg (23.1lb)	11kg (24.2lb)	23kg (57.3lb)	23kg (57.3lb)
Installation	Bench/wall	Bench/wall	Bench/wall	Bench/wall	Bench/wall

Feedwater requirement

Source	Preferably	ally from potable supply, then pre treated. Potable Tap Water hby reverse osmosis (RO) or filtered service deionization (SDI) or distilled.			Originally from potable supply, then pre treated. Preferably reverse osmosis (RO) or filtered service deionization (SDI) or distilled.	
Conductivity	<1 µS/cm ³	<1 µS/cm ³	<1 µS/cm ³	<2000 µS/cm ⁴	<30 µS/cm	
Contaminant						
Hardness	N/A	N/A	N/A	<350 ppm as CaCO ₃	TBC	
Free Chlorine	<0.05 ppm Cl ₂	<0.05 ppm Cl ₂	<0.05 ppm Cl ₂	<0.5 ppm Cl ₂	<0.05 ppm Cl ₂	
Chloramine	<0.05 ppm Cl ₂	<0.05 ppm Cl ₂	<0.05 ppm Cl ₂	<0.2 ppm Cl ₂	<0.05 ppm Cl ₂	
Total Chlorine	<0.05 ppm Cl ₂	<0.05 ppm Cl ₂	<0.05 ppm Cl ₂	<0.5 ppm Cl ₂	<0.05 ppm Cl ₂	
Silica	<2 ppm	<2 ppm	<2 ppm	<30 ppm SiO ₂	<2 ppm SiO ₂	
Carbon Dioxide Maximum	<0.1 ppm	<0.1 ppm	<0.1 ppm	<30 ppm	<30 ppm	
				(recommended <20 ppm)	(recommended <20 ppm)	
Fouling index	<1	<1	<1	<10	<1	
Iron / Manganese	N/A	N/A	N/A	<0.05 ppm	N/A	
Organics (TOC)	N/A	N/A	N/A	<2 ppm	<50 ppb recommended	
Particulates	A 0.2 micron membrane pre filter is recommended for all non-RO feeds to extend point-of-use filter life				N/A	
Temperature	4 - 40°C (Recommended 10 - 25°C)					
Flowrate (requirement at 15°C)	<2 l/min	<2 l/min	<2 l/min	Up to 75 l/hr	Up to 60 l/hr	
Drain requirements (gravity fall with air gap)	None required	None required	None required	>90L/hr	>70L/hr	

³ If feedwater is in the range of 1 µS/cm – 30 µS/cm you will need a pre-conditioning cartridge ⁴ Purification pack life may vary with feedwaters >1400 µS/cm

Feedwater pressure

Maximum	22psi (1.5bar)	22psi (1.5bar)	22psi (1.5bar)	90psi (6bar) ⁵	90psi (6bar)
Minimum	Flooded suction	Flooded suction	Flooded suction	30psi (2bar) ⁵	1psi (0.0 7bar)
516.00 ((1))					

⁵ If <60psi (4 bar) a separate boost pump is recommended

Electrical requirements

Mains Input	100-240V ac, 50-60Hz				
System control voltage (not including pumps and UV)			24V dc		
Power consumption during peak demand	120VA	120VA	120VA	100VA	100VA
Noise level			<40dBa		

ELGA LabWater

Tel: +44 (0) 1494 887500 Fax: +44 (0) 1494 887505 Email: info@elgalabwater.com Website: www.elgalabwater.com

ELGA is the global laboratory water brand name of Veolia Water Solutions and Technologies. The information contained in this document is the property of VWS (UK) Ltd, trading as ELGA LabWater, and is supplied without liability for errors or omissions.

© VWS (UK) Ltd. 2010 – All rights reserved. ELGA® and PURELAB® are registered trademarks of VWS (UK) Ltd.